

ARMY RESEARCH LABORATORY



Fourier Transform Photoacoustic Infrared Spectroscopy of Propellant Formulations, Propellant Ingredients, and Energetic Materials

Jeffrey M. Widder
Kevin L. McNesby

ARL-TR-1232

November 1996

DTIC QUALITY INSPECTION

19961205 021

APPROVED FOR PUBLIC RELEASE; DISTRIBUTION IS UNLIMITED.

NOTICES

Destroy this report when it is no longer needed. DO NOT return it to the originator.

Additional copies of this report may be obtained from the National Technical Information Service, U.S. Department of Commerce, 5285 Port Royal Road, Springfield, VA 22161.

The findings of this report are not to be construed as an official Department of the Army position, unless so designated by other authorized documents.

The use of trade names or manufacturers' names in this report does not constitute indorsement of any commercial product.

REPORT DOCUMENTATION PAGE			Form Approved OMB No. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project(0704-0188), Washington, DC 20503.				
1. AGENCY USE ONLY (Leave blank)		2. REPORT DATE November 1996		3. REPORT TYPE AND DATES COVERED Final, Jun - Dec 95
4. TITLE AND SUBTITLE Fourier Transform Photoacoustic Infrared Spectroscopy of Propellant Formulations, Propellant Ingredients, and Energetic Materials			5. FUNDING NUMBERS PR: 1L162618AH80	
6. AUTHOR(S) Jeffrey M. Widder and Kevin L. McNesby				
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) U.S. Army Research Laboratory ATTN: AMSRL-WT-PC Aberdeen Proving Ground, MD 21005-5066			8. PERFORMING ORGANIZATION REPORT NUMBER ARL-TR-1232	
9. SPONSORING/MONITORING AGENCY NAMES(S) AND ADDRESS(ES)			10. SPONSORING/MONITORING AGENCY REPORT NUMBER	
11. SUPPLEMENTARY NOTES				
12a. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release; distribution is unlimited.			12b. DISTRIBUTION CODE	
13. ABSTRACT (Maximum 200 words) Infrared spectra of 48 energetic materials, propellant ingredients, and propellant formulations have been measured using a Fourier transform infrared (FTIR) spectrometer employing photoacoustic (PA) detection. Spectra are divided among the categories of propellant formulations, neat energetic materials, binders, plasticizers, stabilizers and others. Spectra of propellant formulations have been annotated, where appropriate, to indicate features due to individual ingredients.				
14. SUBJECT TERMS spectra, propellant formulations, energetic materials, photoacoustic spectroscopy			15. NUMBER OF PAGES 67	
			16. PRICE CODE	
17. SECURITY CLASSIFICATION OF REPORT UNCLASSIFIED	18. SECURITY CLASSIFICATION OF THIS PAGE UNCLASSIFIED	19. SECURITY CLASSIFICATION OF ABSTRACT UNCLASSIFIED	20. LIMITATION OF ABSTRACT UL	

INTENTIONALLY LEFT BLANK.

PREFACE

This work was performed while the author held a National Research Council-Army Research Laboratory Research Associateship.

INTENTIONALLY LEFT BLANK.

TABLE OF CONTENTS

		<u>Page</u>
	PREFACE	iii
1.	INTRODUCTION	1
2.	BACKGROUND	1
3.	EXPERIMENTAL	3
4.	DISCUSSION AND RESULTS	3
5.	CONCLUSIONS	5
6.	REFERENCES	7
	APPENDIX A: SPECTRA OF FORMULATIONS	9
	APPENDIX B: SPECTRA OF NEAT ENERGETIC MATERIALS	23
	APPENDIX C: SPECTRA OF BINDERS	33
	APPENDIX D: SPECTRA OF PLASTICIZERS	45
	APPENDIX E: SPECTRA OF STABILIZERS	53
	APPENDIX F: SPECTRA OF OTHERS	59
	LIST OF ABBREVIATIONS	69
	DISTRIBUTION LIST	71

INTENTIONALLY LEFT BLANK.

1. INTRODUCTION

Propellant formulations and their ingredients are often opaque and/or poorly reflective, making it difficult to measure infrared (IR) spectra using conventional techniques. Techniques for measuring IR spectra, which rely on detection of beam attenuation, such as transmission, reflectance, and attenuated total reflectance (ATR), are often not suitable for such samples or require elaborate sample preparation. To avoid saturation, transmission spectroscopy in the mid-IR spectral region requires that the path length through a typical propellant sample be on the order of 5–10 μm . Reflectance spectroscopy of many crystalline energetic materials may be sensitive to sample alignment and frequently requires mathematical corrections to the raw data. ATR spectroscopy requires that the sample be brought into intimate contact with an internal reflection element and may not be suitable for many solid samples. Fourier transform infrared (FTIR) spectroscopy employing photoacoustic (PA) detection may be used to overcome many of the limitations posed by these traditional techniques. PA detection requires minimal sample preparation and may be successfully used with opaque, poorly reflecting samples and samples with widely varying morphology (McClelland et al. 1993). Additionally, for some samples, the PA technique allows measurement of vibrational spectra at discrete depths beneath the sample surface (Palmer 1993).

This report contains a PA IR spectral library of propellant formulations, neat ingredients, and energetic materials. At the present time, we know of no commercially available IR libraries dealing specifically with spectra of propellant formulations, propellant ingredients, and energetic materials. We believe such a library would have applications in and beyond those of the propellant community, especially in the areas of environmental site remediation, identification of unknowns, and the analysis of stored propellants that may be in varying states of decomposition.

2. BACKGROUND

Following the excitation of vibrational modes by the absorption of IR radiation, relaxation processes dissipate the absorbed energy. Most of the energy is distributed among the normal modes according to the Boltzmann distribution, resulting in a temperature rise of the sample. If the intensity of the incident IR radiation is constant, the rate of heat generation in the sample will be constant. If, however, the intensity of the incident IR radiation is modulated, the heating of the sample will be modulated. The modulation of the sample surface temperature causes the gas adjacent to the sample surface to expand at the modulation frequency of the incident radiation. The modulated expansion of the gas adjacent to the

sample causes pressure pulses in the gas. When the modulation frequency is in the acoustic range, the pressure pulses may be detected using a microphone. The magnitude of the signal reaching the microphone, during one modulation period, is proportional to the strength of the sample absorption. When the modulation of the IR radiation is accomplished using an interferometer, a Fourier transform of the microphone output yields the IR spectrum of the sample (Dittmar, Palmer, and Carter 1994).

PA signal generation is dependent on the transfer of vibrational energy from the sample to translational energy of the gas adjacent to the sample surface. This results in a time lag (phase shift) between IR absorption by the sample and detection of signal by the microphone. The magnitude of the time lag is dependent on the depth within the sample at which absorption occurs. This time lag between when the IR radiation is absorbed and when the resulting sound is detected by the microphone is manifested as a phase difference between IR source modulation frequency and microphone output frequency. The significance of the phase difference is that it can be used to calculate the depth beneath the sample surface of the absorbing species (Palmer 1993). At high modulation frequencies, only heat produced by absorption on or near the surface contributes to the signal detected during the modulation period. At lower modulation frequencies, the detected signal contains contributions from absorption deeper within the sample. Comparison of spectra measured at different modulation frequencies can reveal concentration profiles and/or stratification of absorbing species in the sample surface region (Dittmar 1994).

IR spectral depth profiling by PA detection is most easily accomplished using a step-scan interferometer employing phase-sensitive lock-in detection. This is because in the step-scan mode of operation, the modulation of the IR intensity is not wavelength-dependent. The technique involves positioning the moving mirror of the interferometer at a known retardation (optical path difference) and modulating the throughput of the interferometer at this mirror position by either mechanically chopping the beam or by applying a small amplitude sinusoidal "dither" to the moving mirror. The frequency of the modulation (the chopping or dithering frequency) at the nominally fixed-mirror position is then used as the reference frequency by the lock-in amplifier. In rapid scan interferometry, the IR modulation frequency and, thus, the effective sampling depth at which absorption can be detected is wavelength-dependent. The ability to preserve spectral multiplexing while removing the optical retardation dependence of the wavelength modulation, achievable by step-scan operation (Palmer 1993), allows uniform depth profiling over the full spectral range of the instrument. By varying the phase at which signal is collected relative to the chopping or dithering frequency, a single spectrum may be employed to yield absorption spectra at several discrete depth intervals beneath the sample surface (Palmer 1993).

3. EXPERIMENTAL

All the spectra in this library were measured using an MTEC 100 PA accessory and a Mattson Polaris spectrometer operating at a spectral resolution of 4 cm^{-1} . The spectra are the result of 32 or 64 coadded scans of the interferometer and have been ratioed to a spectrum of carbon black powder measured on the same day. Ratioing the spectra to the spectrum of carbon black normalizes the spectra to the instrumental throughput and the detector response (Palmer 1993). This is equivalent to ratioing a transmission spectrum to the spectrum of the open beam. Samples were placed in the sample holder of the MTEC accessory (a 6-mm-diameter cylindrical cup approximately 3 mm in depth), and the accessory was purged with dry helium gas to remove water vapor and CO_2 gas. An additional benefit of purging with helium is that the low heat capacity of the gas gives a signal intensity two to three times that of a cell filled with dry air.

Samples of propellant formulations, in the form of thin cross sections of grains and cords, were prepared by cutting with a razor blade. Samples of propellant formulations, when in the form of sheets, were prepared using a small sheet-metal punch. Samples of neat ingredients and energetic materials, when in powder form, were added to the sample cup using a small spatula and metal funnel. Samples of neat ingredients and energetic materials, when in liquid form, were added to the sample cup with a disposable pipet.

4. DISCUSSION AND RESULTS

The spectral library is divided into six categories. These categories are propellant formulations (Appendix A), energetic materials (Appendix B), binders (Appendix C), plasticizers (Appendix D), stabilizers (Appendix E), and others (Appendix F). The category of propellant formulations has been further broken down into four subcategories based on the formulation type. These formulation types are single-base propellants (nitrocellulose-based), double-base propellants (nitrocellulose- and nitroglycerin-based), triple-base propellants (nitrocellulose-, nitroglycerin-, and nitroguanidine-based), and nitramine-based (usually RDX- or HMX-based, or a combination of double base and nitramine-based). A list of the propellant formulations for which photoacoustic Fourier-transform infrared (PA-FTIR) spectra are reported here, arranged by propellant type, is given in Table 1. Table 2 lists the neat ingredients for which spectra are reported. Table 3 shows the major components of the propellant formulations listed in Table 1. Note that the List of Abbreviations section at the end of the report provides all the chemical abbreviations along with their spelled-out versions.

Table 1. Propellant Formulations

Single Base	Double Base	Triple Base	Nitramine-Based
M10	JA2 RPD351 M44, M9	M30	C4, XM39, M43, JAG, JAX(2R20), HELP19, JAX3

Table 2. Neat Ingredients

Energetic	Binder	Plasticizer	Stabilizer	Other
NC	NC	ATEC	AKARDIT II	KB ₁₁ H ₁₄
NG	CAB	Paraplex G-59	EC	TNB
RDX	PDNPA	DEP	MC	DNP
RDX (Chinese)	hycar 4051	DEGDN		DMNP
CL 20	TPE	NG		triazole
HMX	CAB 500-5	triacetin		CNNQ
TNT	hycar 4054			CP
TNAZ	BEMO-BMMO/THF			
NQ	BAMO-BAMO/AMMO-BAMO			

Table 3. Propellant Composition

Propellant	Composition (approximate percentages)
JA2	59.9% NC, 14.9% NG, 24.8% DEGDN, 0.1% AKARDIT II
RPD351	49.5% NC, 40.9% NG, 7.2% DEGDN, 0.8% AKARDIT II, 1.5% KN
JAG	51.2% NC, 18% NG, 11.5% DEGDN, 0.7% AKARDIT II, 3.5% paraplex G-59, 15% RDX
JAX (2R20)	50.6% NC, 11.4% NG, 19.2% DEGDN, 0.7% AKARDIT II, 17.8% RDX
M44	52.9% NC, 43% NG, 2% EC, 2% DEP
M30	28% NC, 22.5% NG, 47.7% NQ, 1.5% EC
M9	57.6% NC, 39.9% NG, 0.75% EC, 1.5% KN, 0.35% ALC
M43	76% RDX, 12% CAB, 7.6% plasticizer, 4% NC, 0.4% EC
XM39	76% RDX, 12% CAB, 7.6% ATEC, 4% NC, 0.4% EC
M10	96% NC, 1% DPA, 1% KS, 1.5% water
HELP19	76% RDX, 11% CAB, 4% DANPE, 4% NC, 0.4% EC

5. CONCLUSIONS

PA-FTIR spectroscopy is an excellent technique for measuring spectra of propellant formulations and neat ingredients of these formulations. The spectra presented in this library are the initial spectra that have been measured for a larger library. These spectra have also been used to train a neural network that will be used to identify unknown formulations and quantify component concentrations for both quality control and aging studies. The result of preliminary training of the neural network were presented at the 1995 Federation of Analytical Chemists and Spectroscopist Societies meeting. This spectral library should be useful for identification of unknowns and also useful as a guide for determination of changes in propellant composition during storage.

A planned conversion of our spectrometer to step-scan operation will allow for spectral depth profiling of propellant surfaces. This will permit measurement of ingredient concentration profiles near the surface of propellant grains, providing valuable information that may be characteristic of manufacture and/or storage condition.

INTENTIONALLY LEFT BLANK.


6. REFERENCES


- Dittmar, R. M., R. A. Palmer, and R. O. Carter III. "Fourier Transform Photoacoustic Spectroscopy of Polymer." Applied Spectroscopy Reviews, vol. 29, p. 171, 1994.
- McClelland, J. F., R. W. Jones, S. Luo, and L. M. Seaverson. "A Practical Guide to FT-IR Photoacoustic Spectroscopy." Practical Sampling Technique for Infrared Analysis, edited by Patricia B. Coleman, ch. 5, 1993.
- Palmer, Richard A. "Photoacoustic and Photothermal Spectroscopies." Determination of Electronic and Optical Properties, edited by Bryant W. Rossiter and Roger C. Baetzold, Physical Methods of Chemistry Series, Second edition, vol. 8, ch. 2, 1993.
- Palmer, Richard A. "Step-Scan FT-IR: A Versatile Tool for Time- and Phase-Resolved Vibrational Spectroscopy." Spectroscopy, vol. 8, no. 2, p. 26, 1993.


INTENTIONALLY LEFT BLANK.


APPENDIX A:
SPECTRA OF FORMULATIONS


INTENTIONALLY LEFT BLANK.

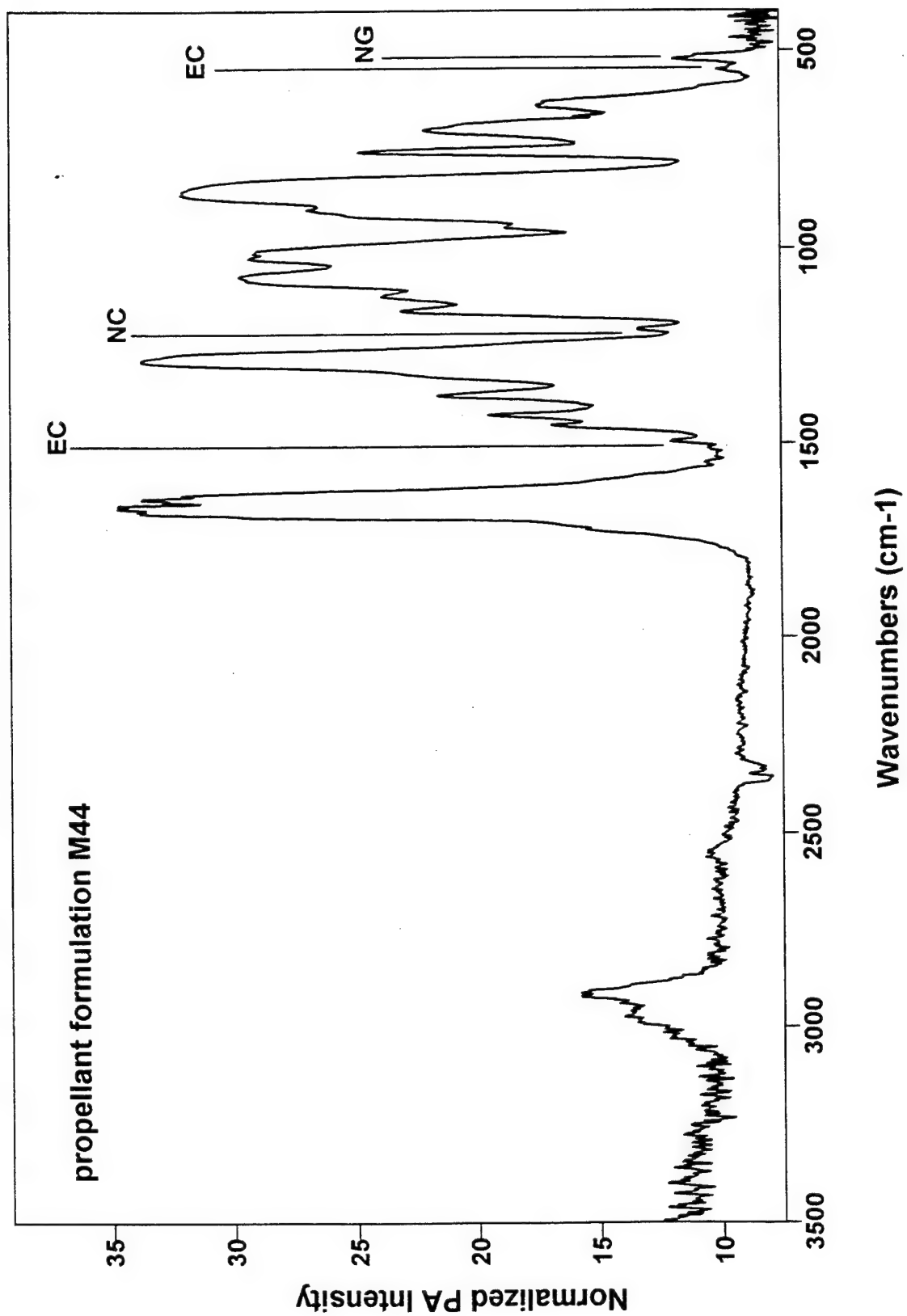


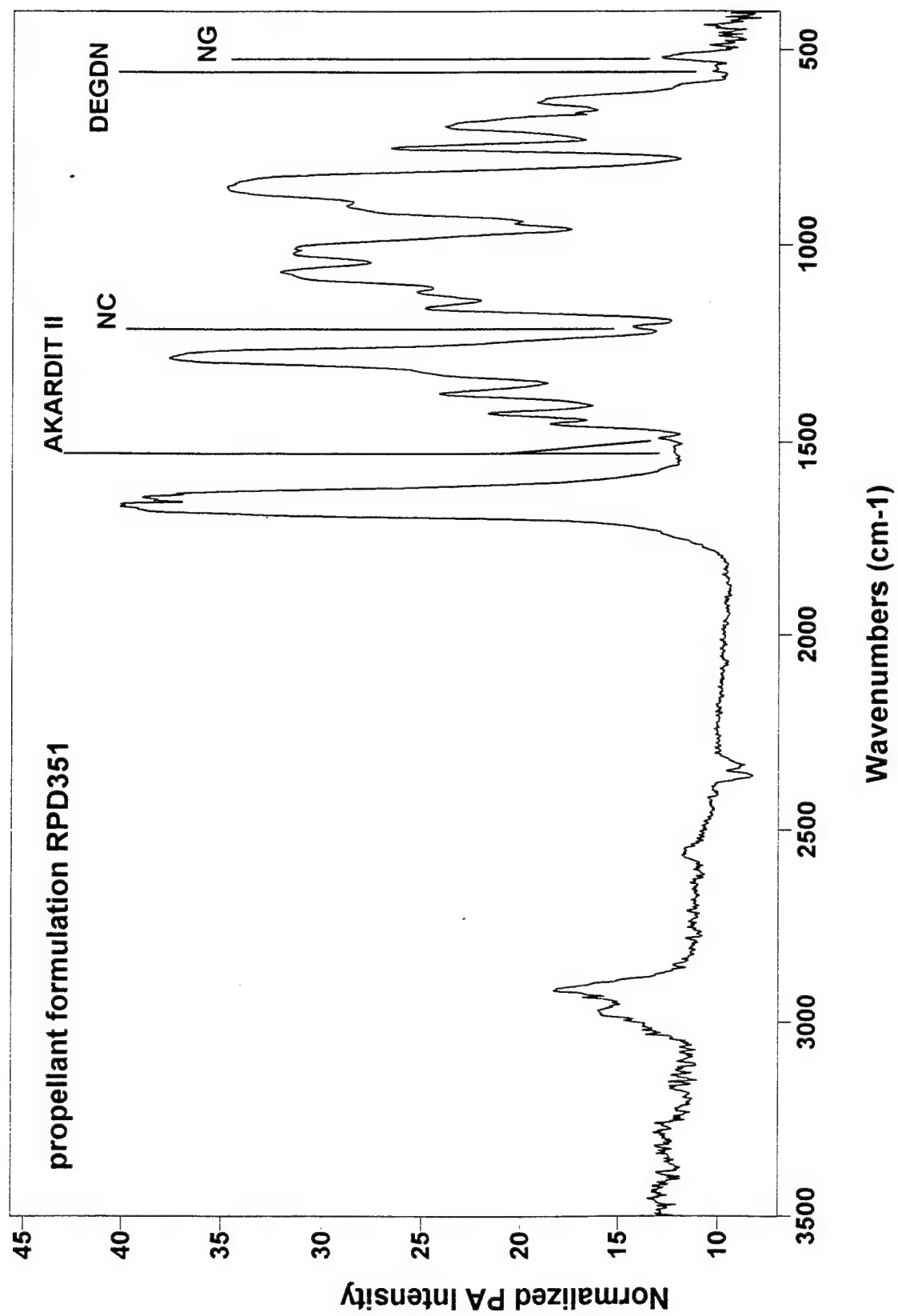


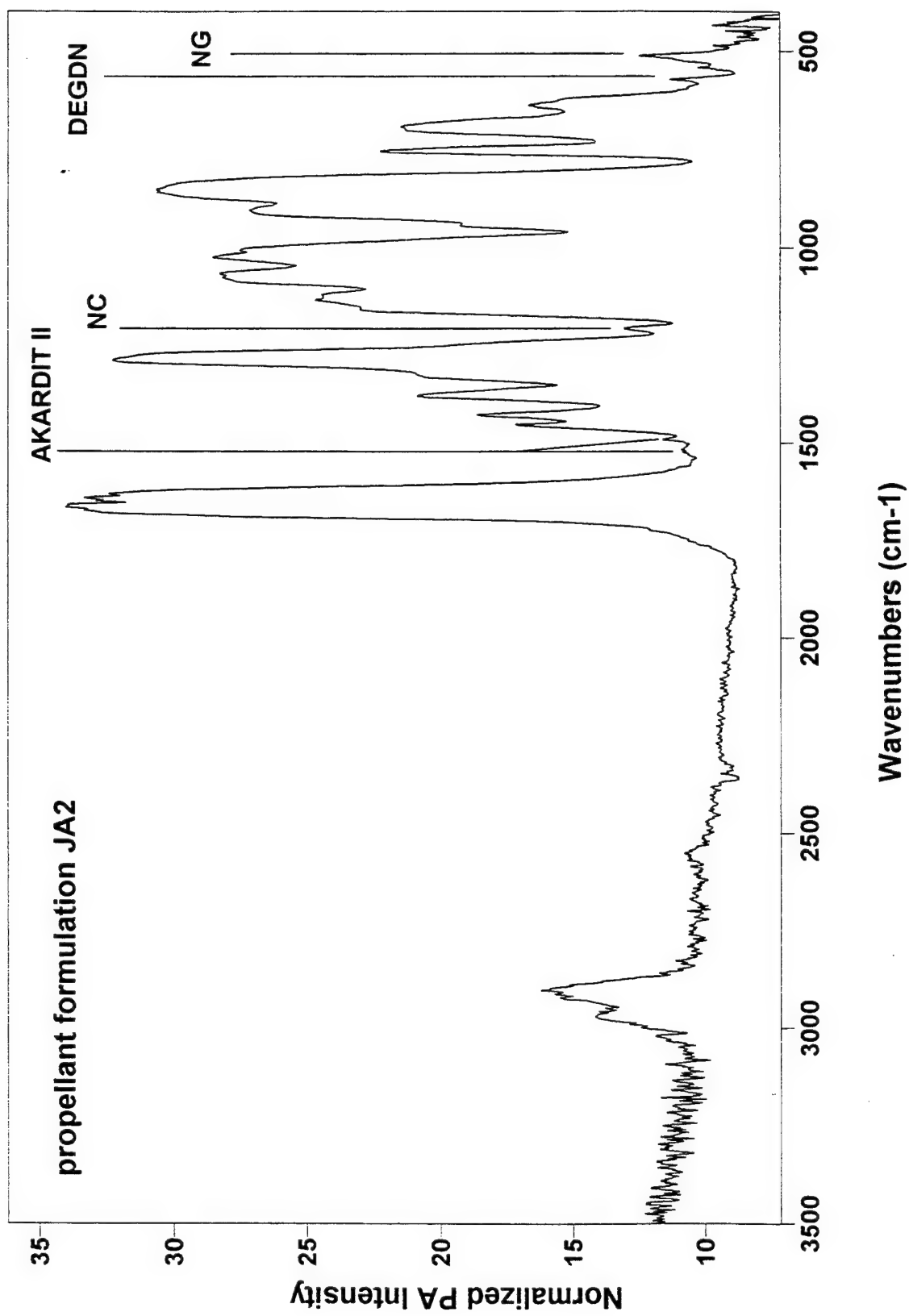


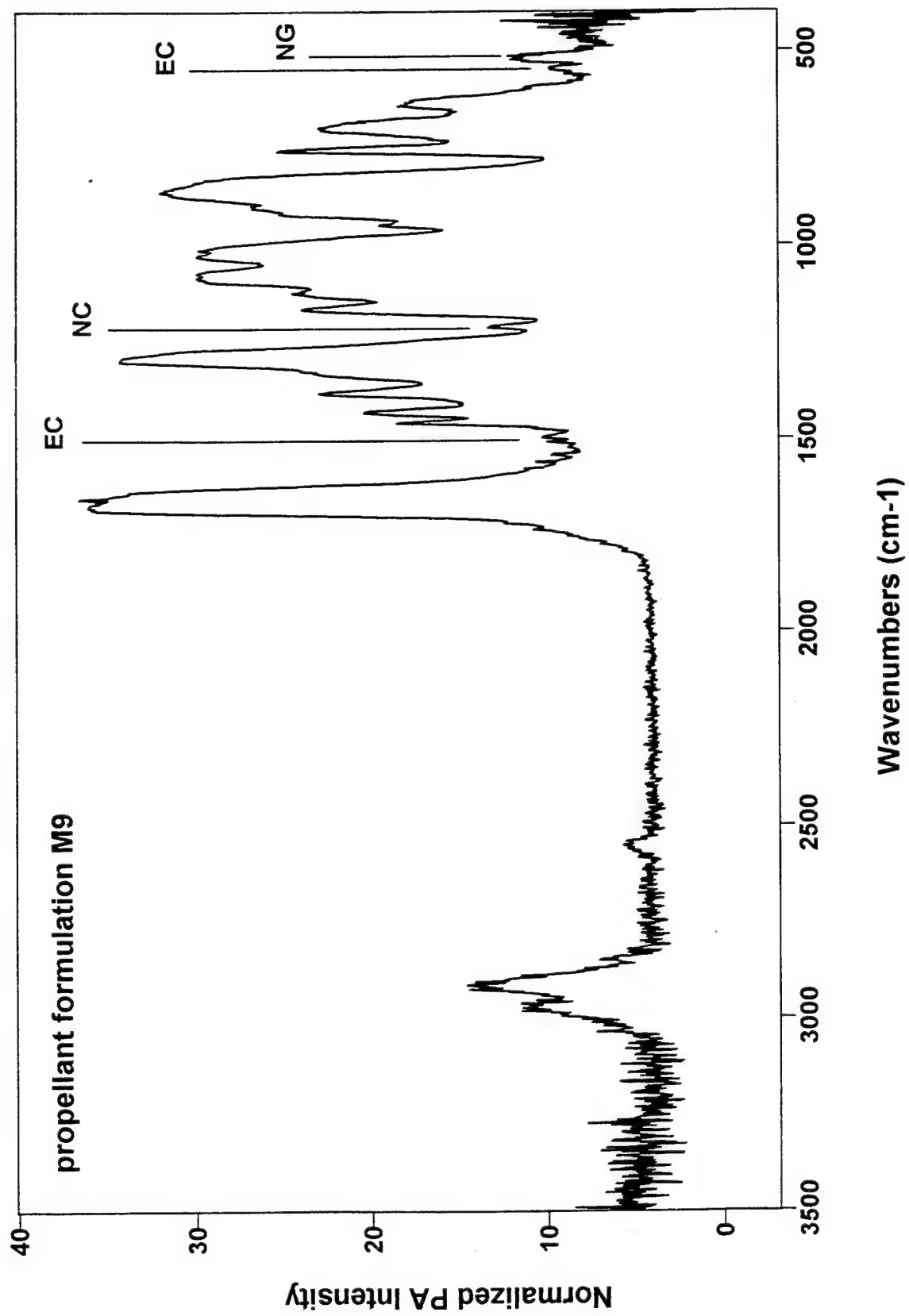


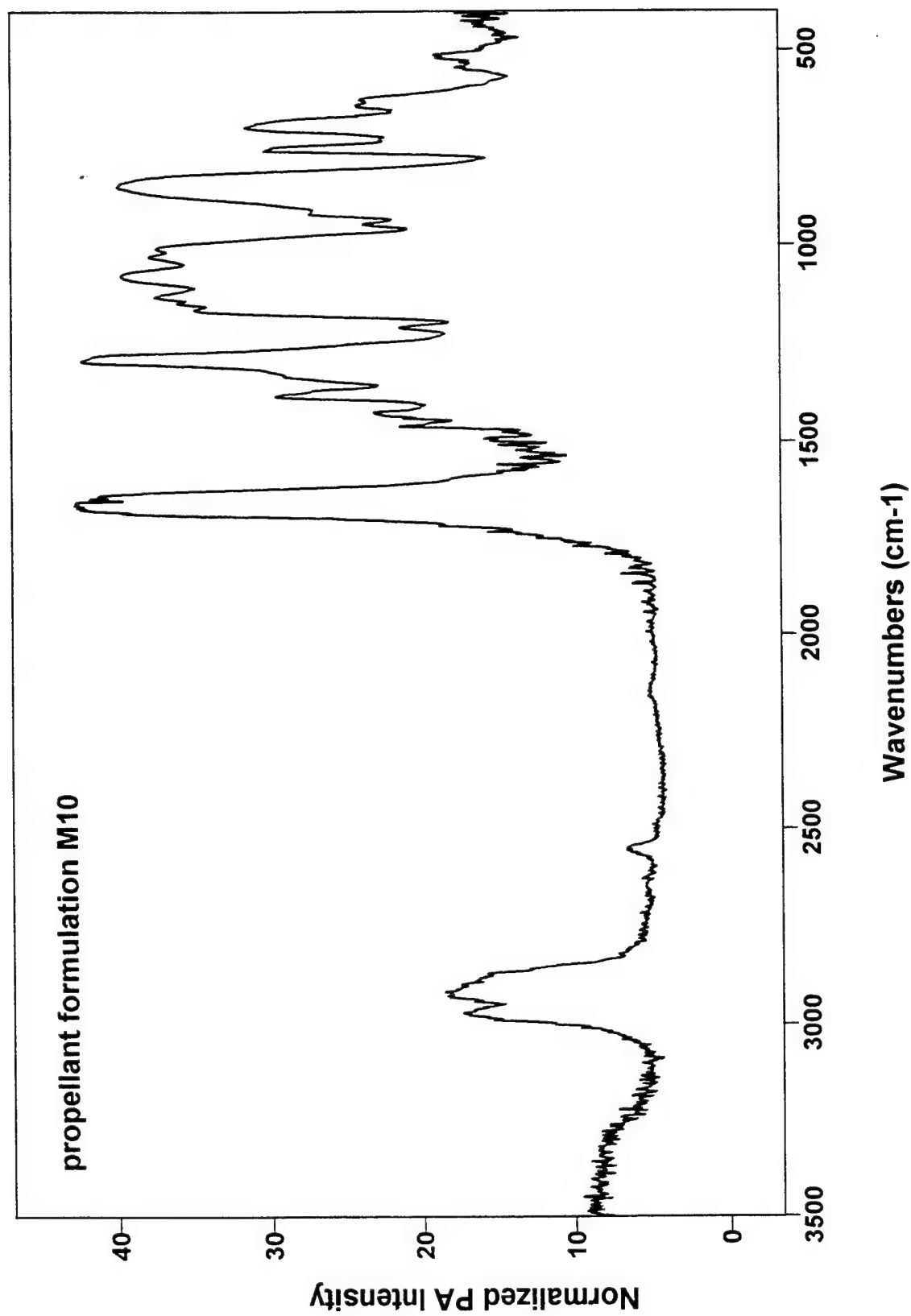


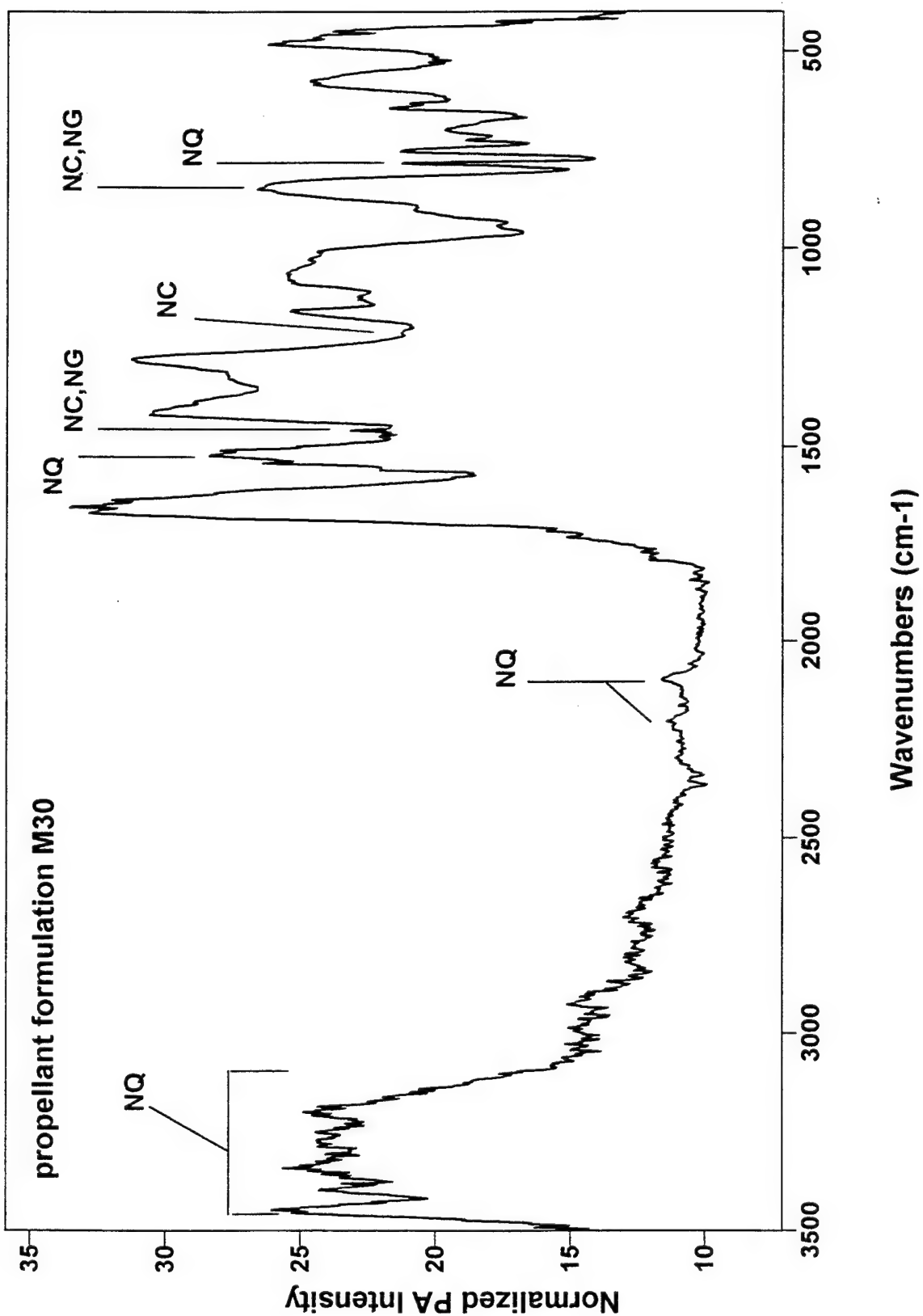


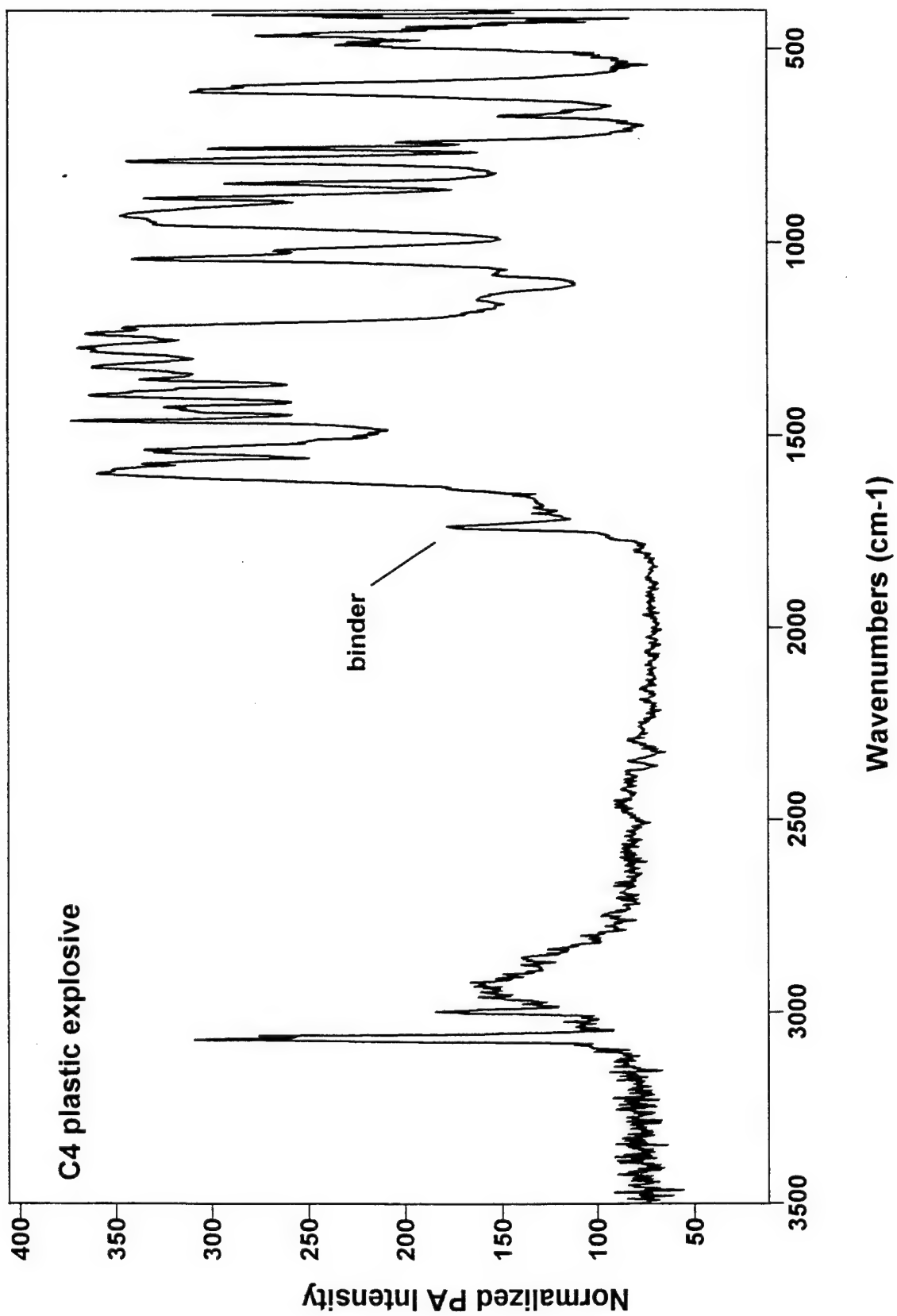


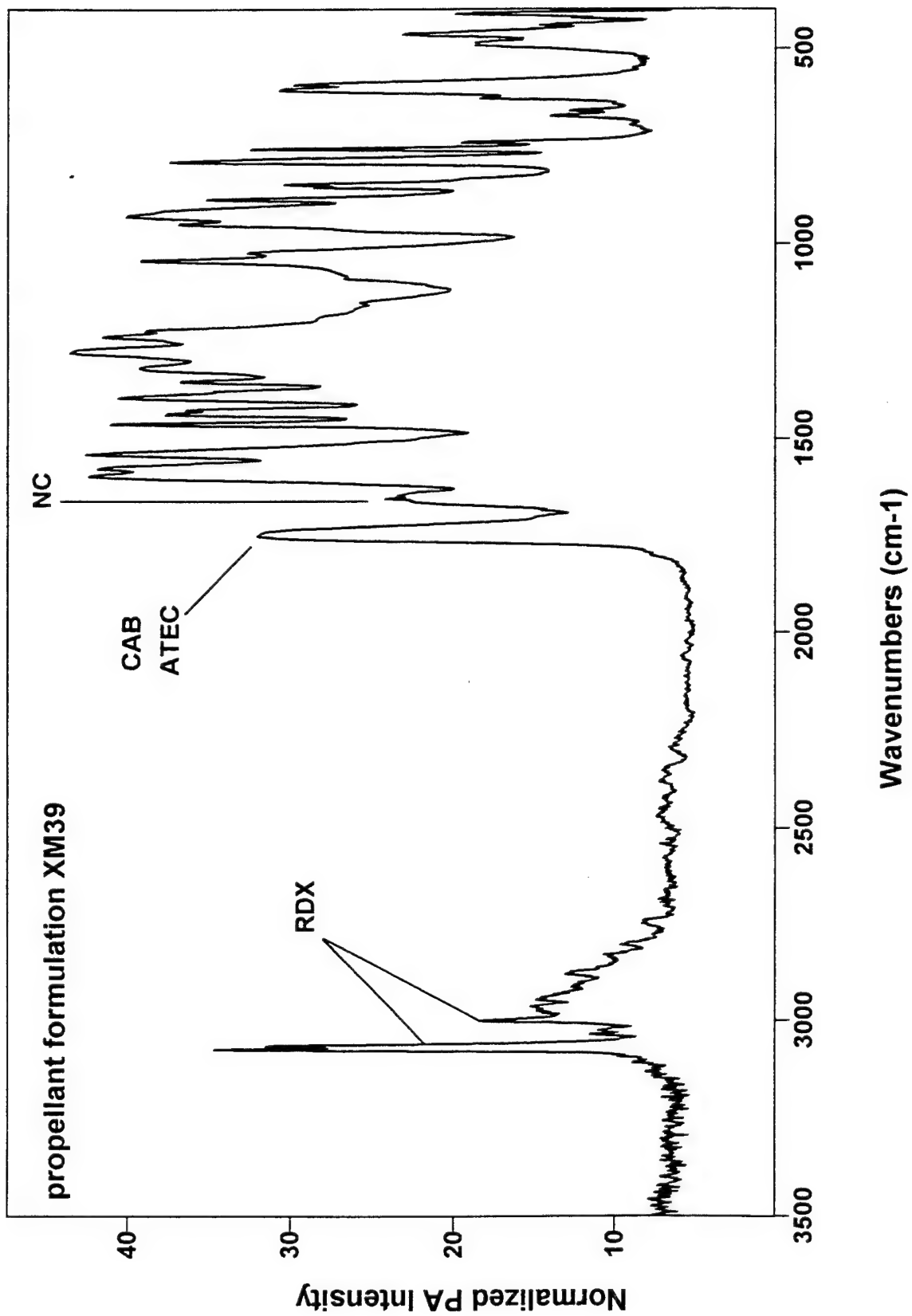

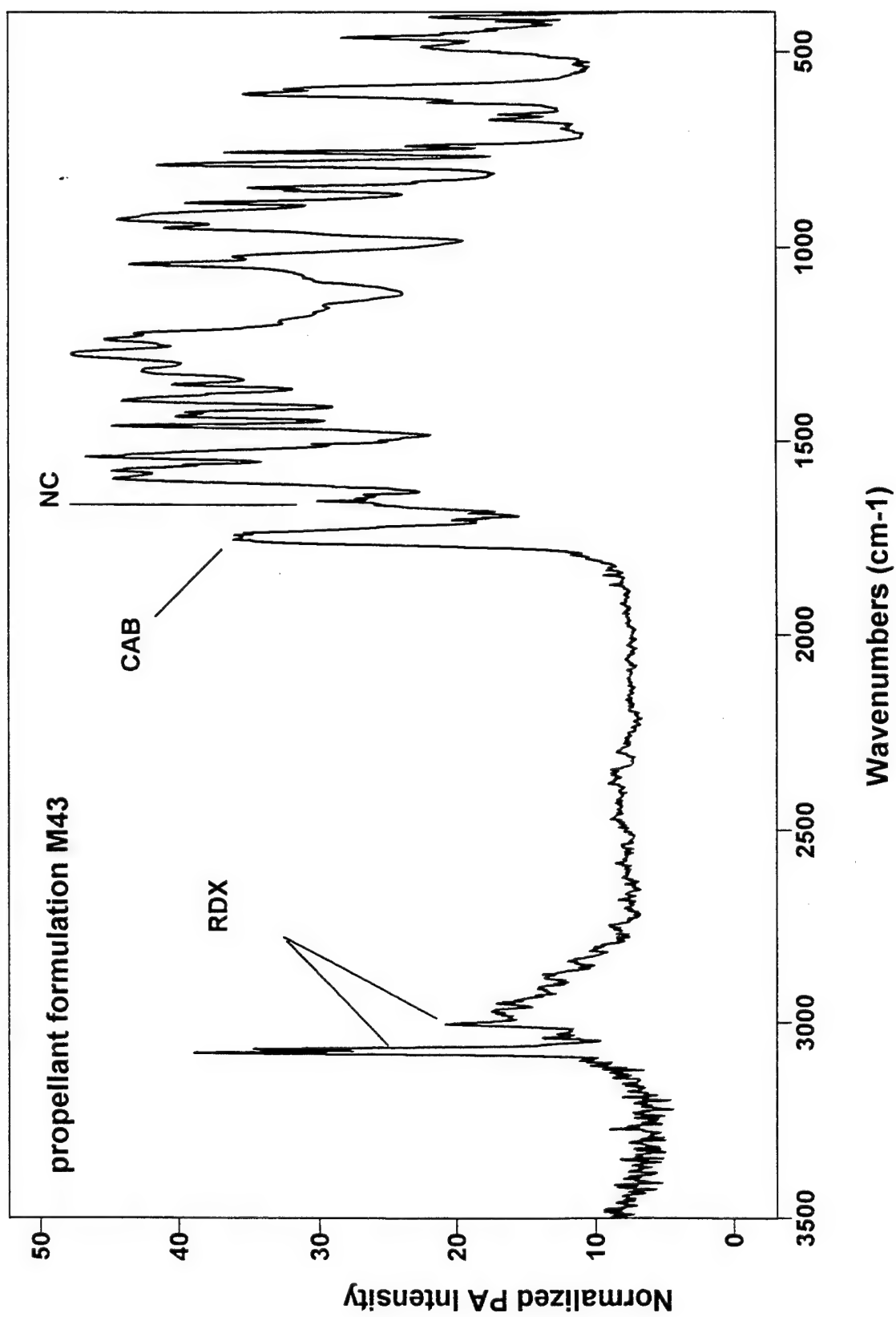

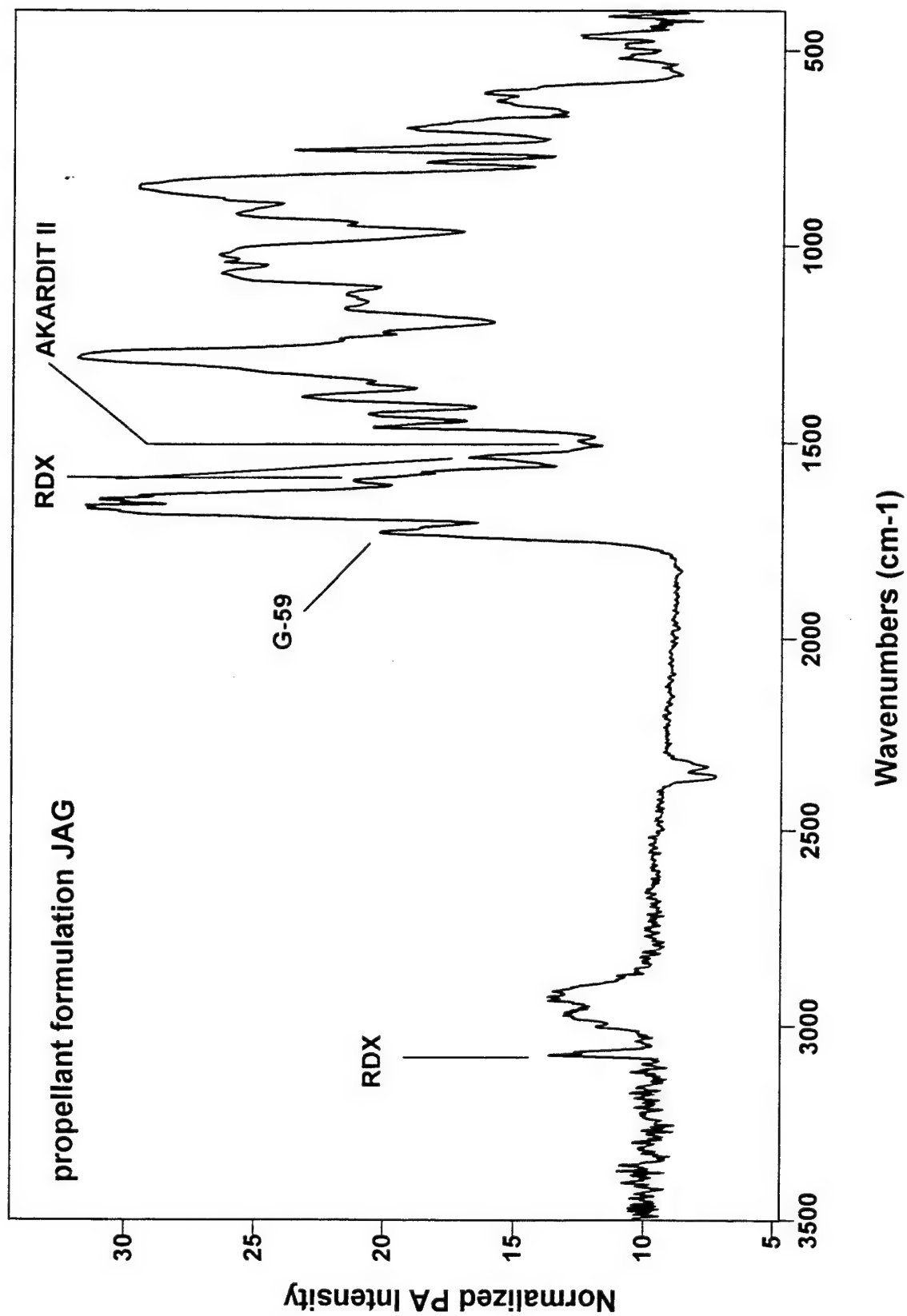


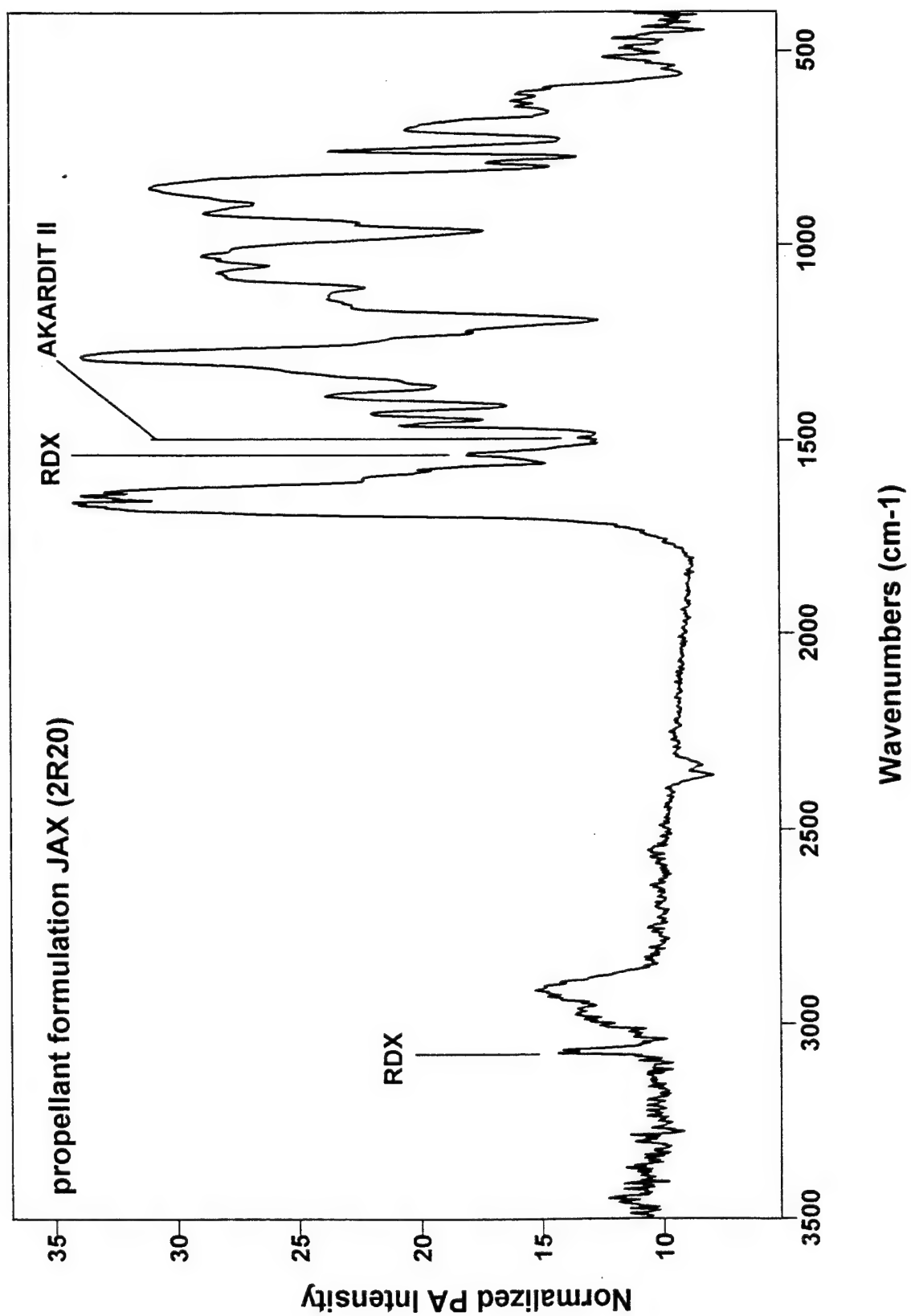


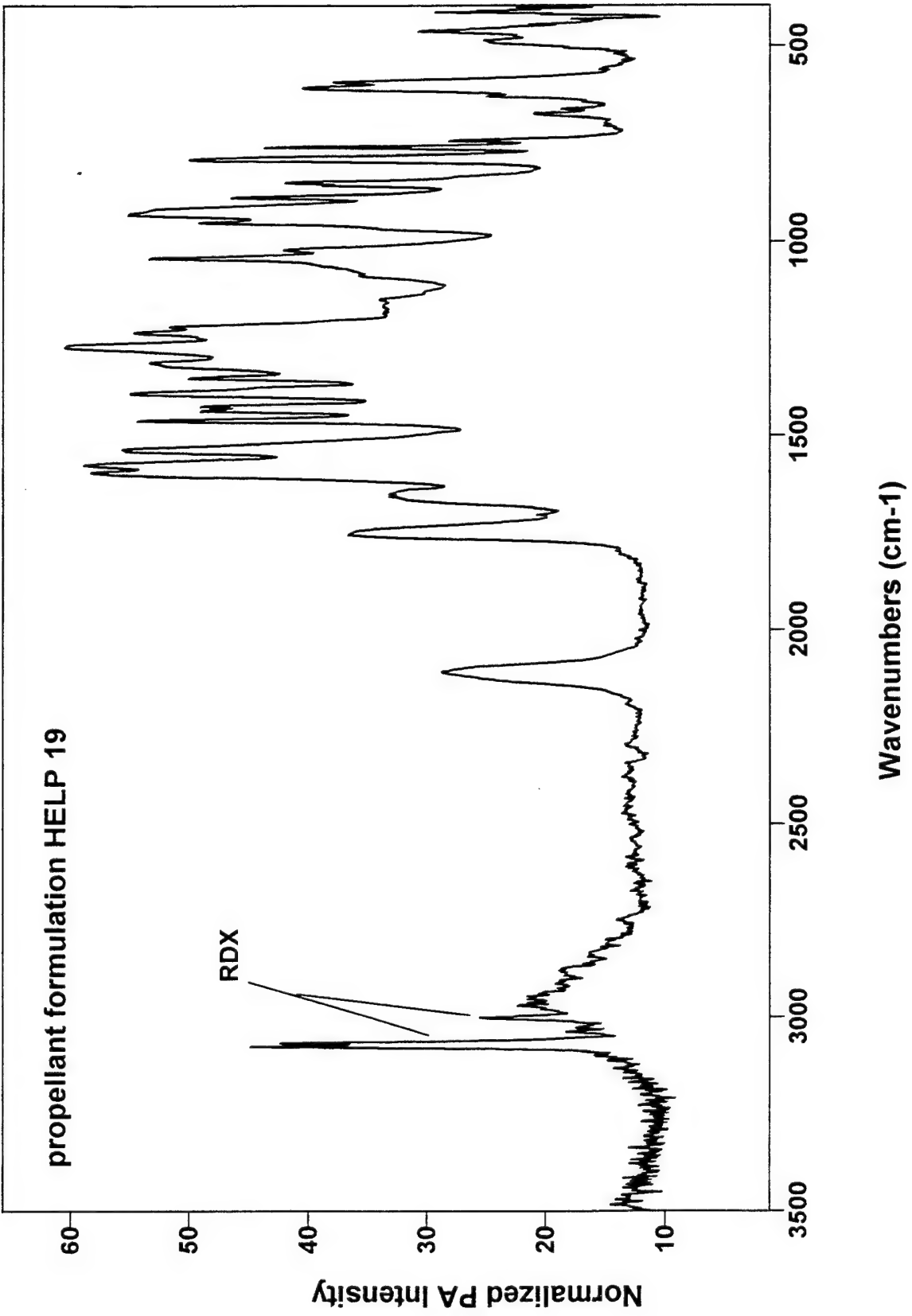



APPENDIX B:
SPECTRA OF NEAT ENERGETIC MATERIALS


INTENTIONALLY LEFT BLANK.

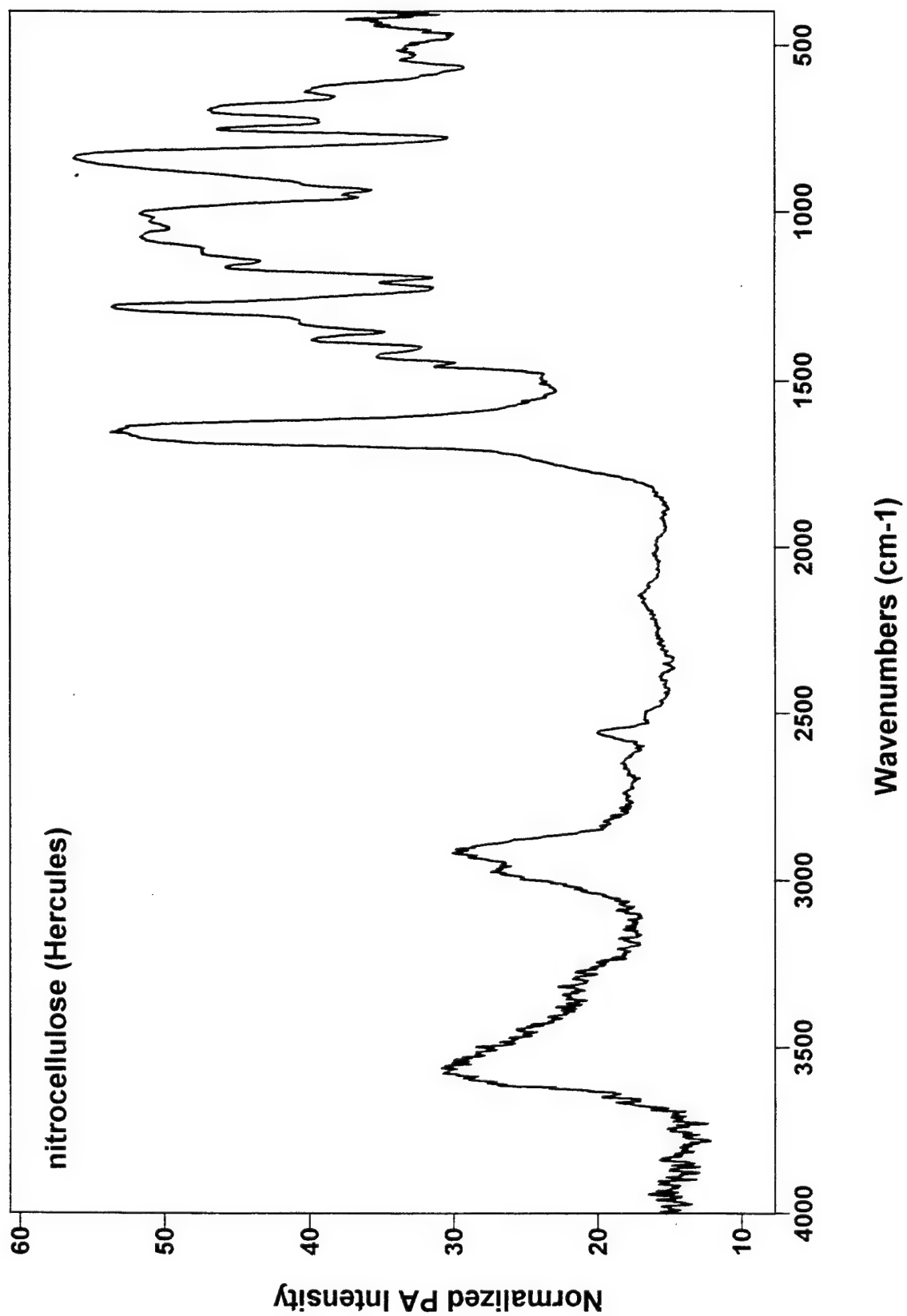


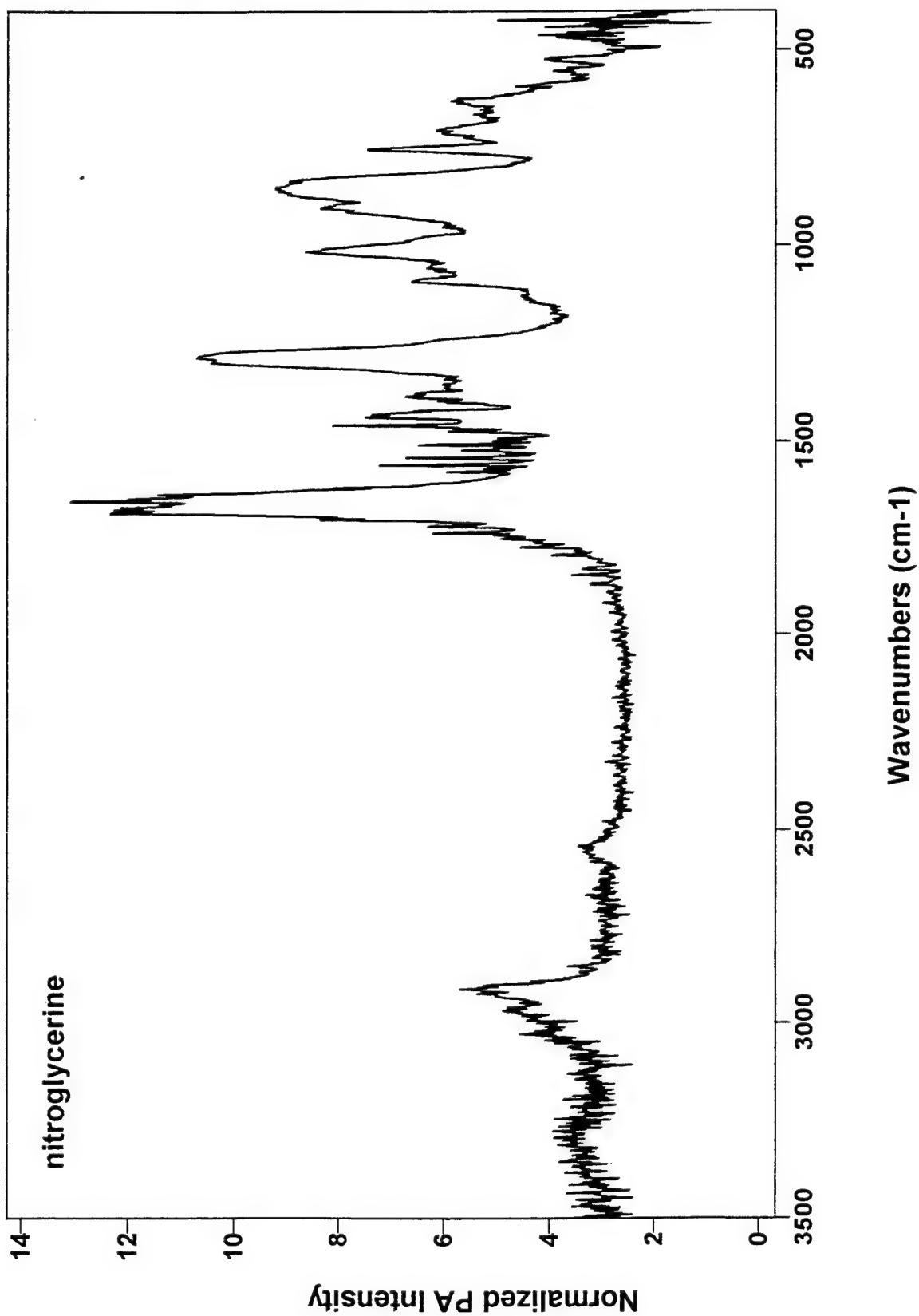


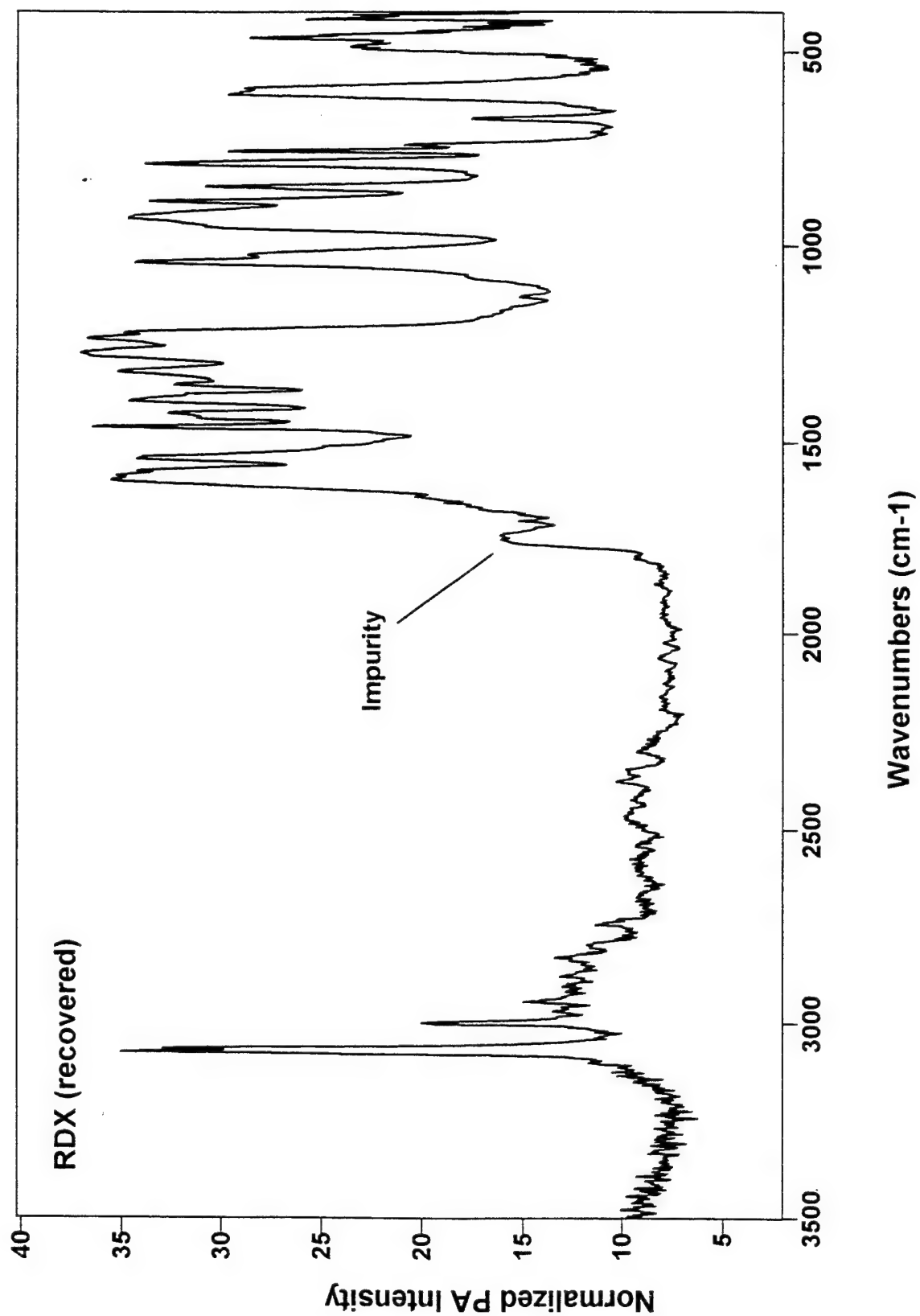


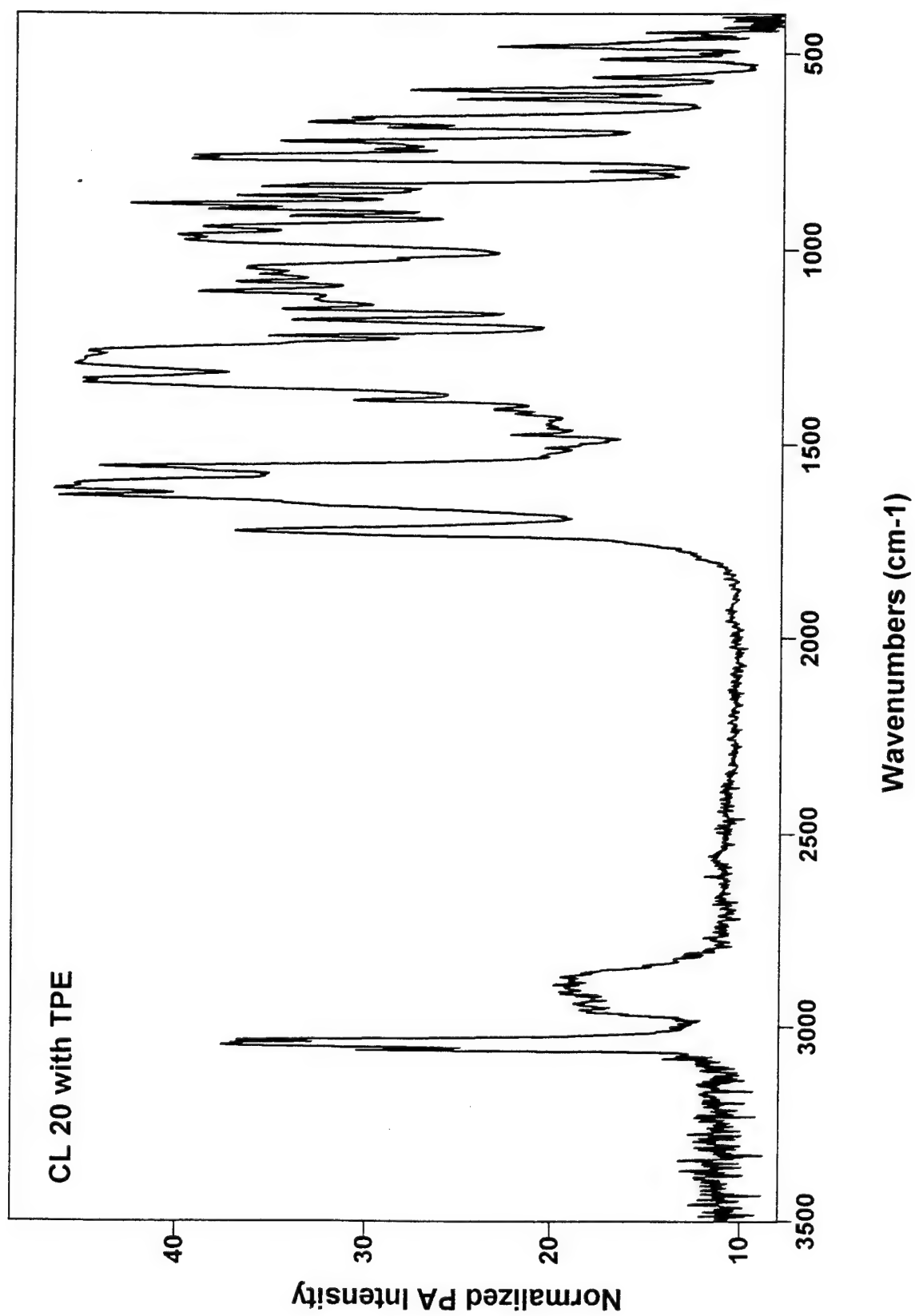

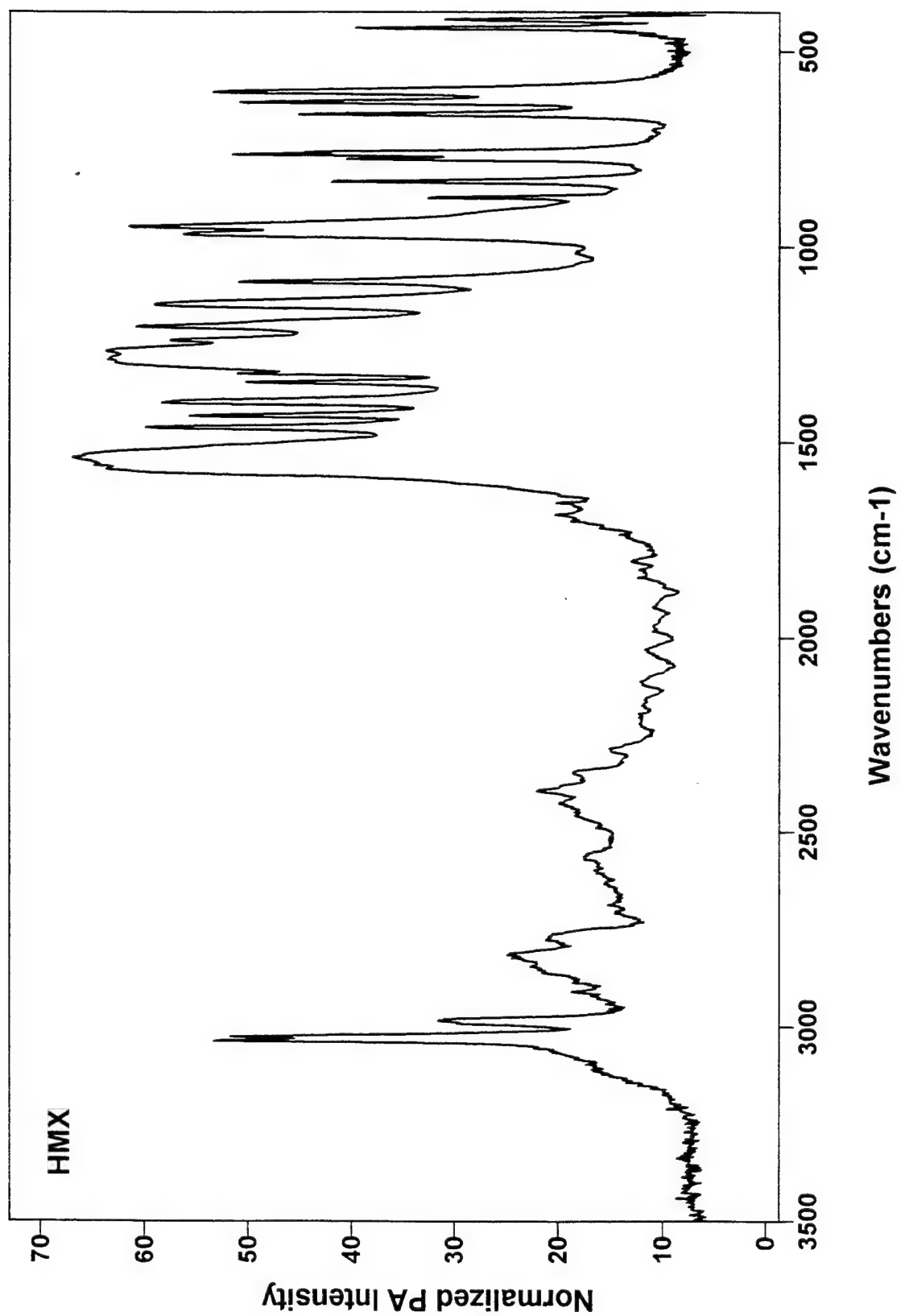

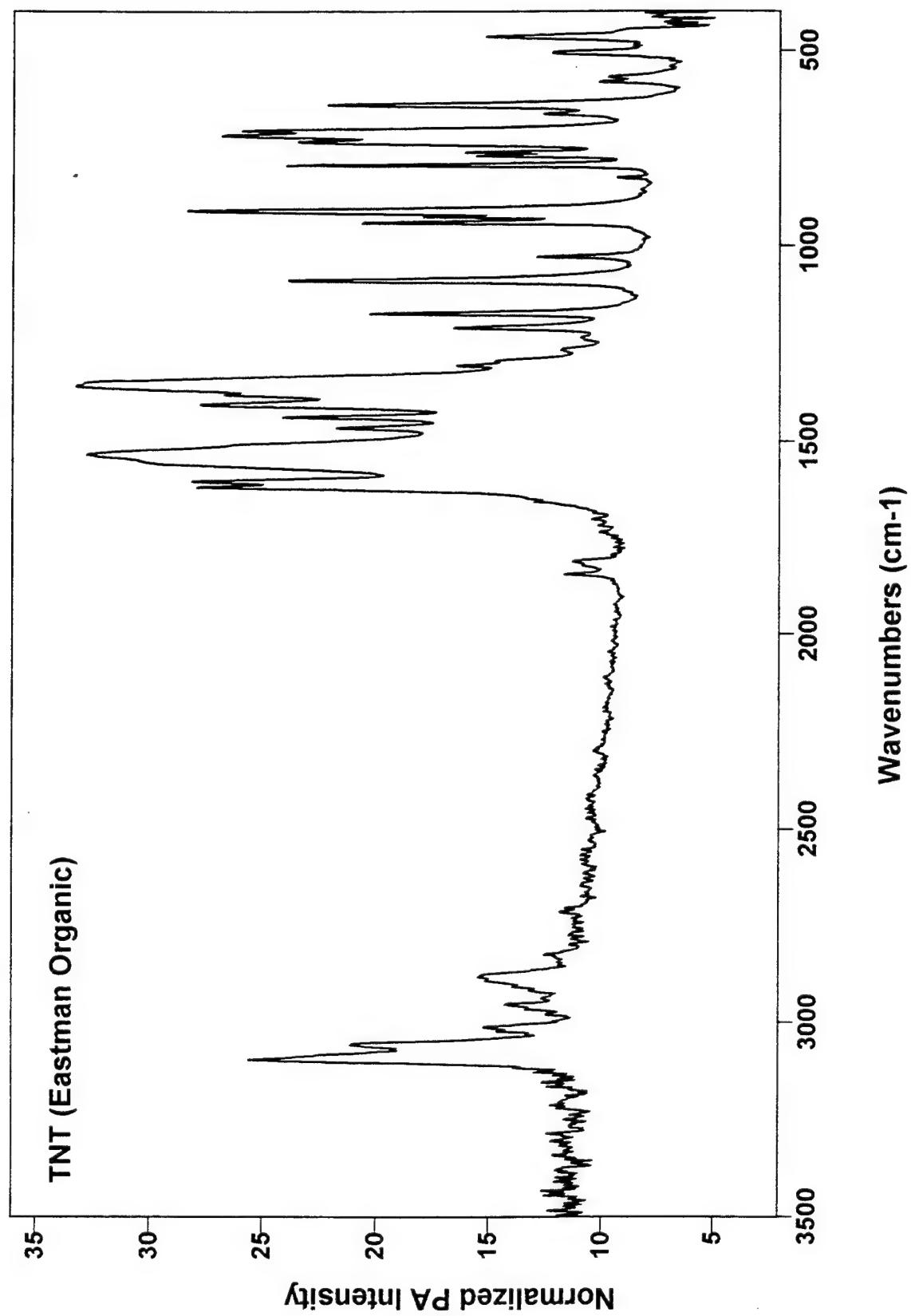


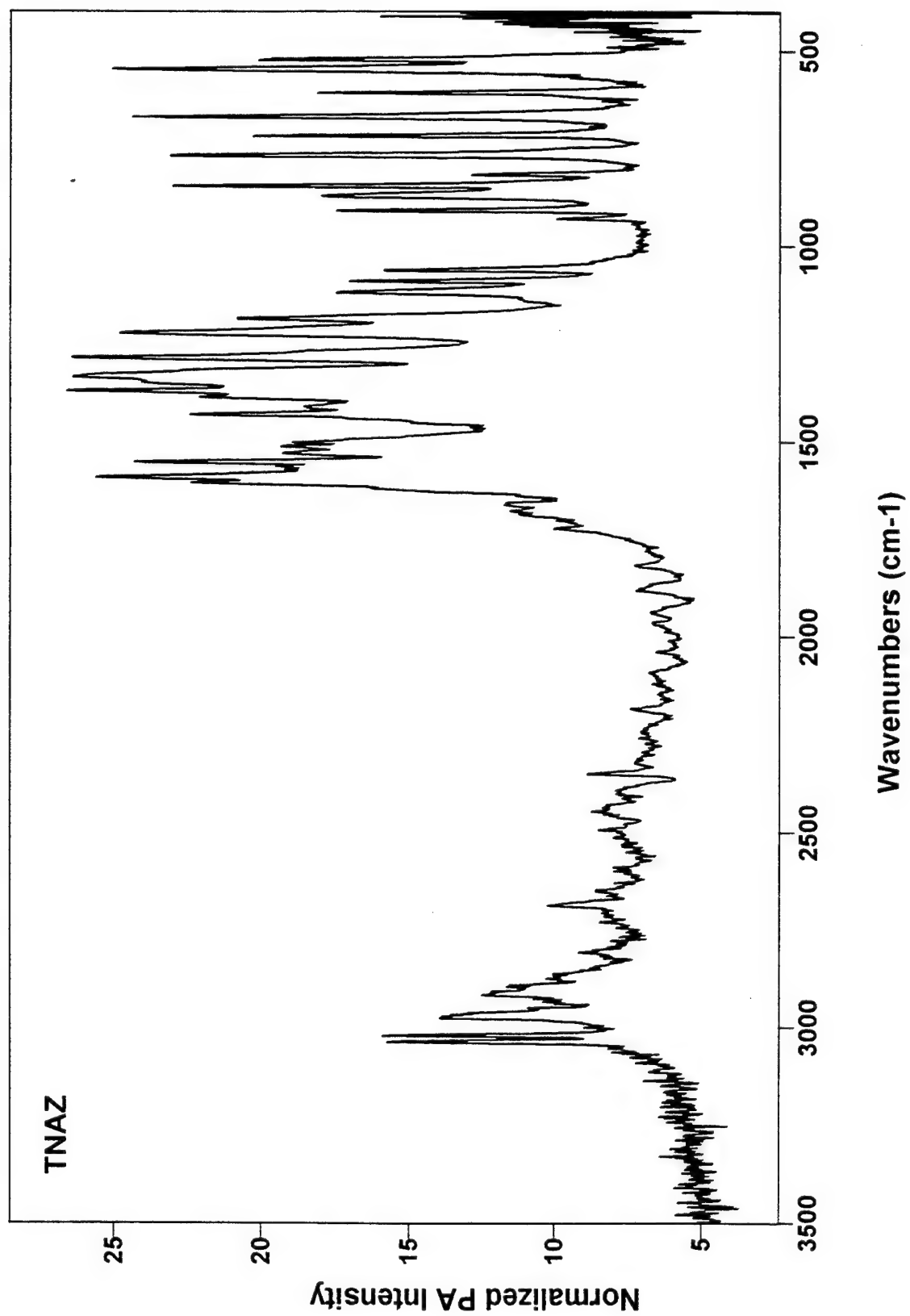


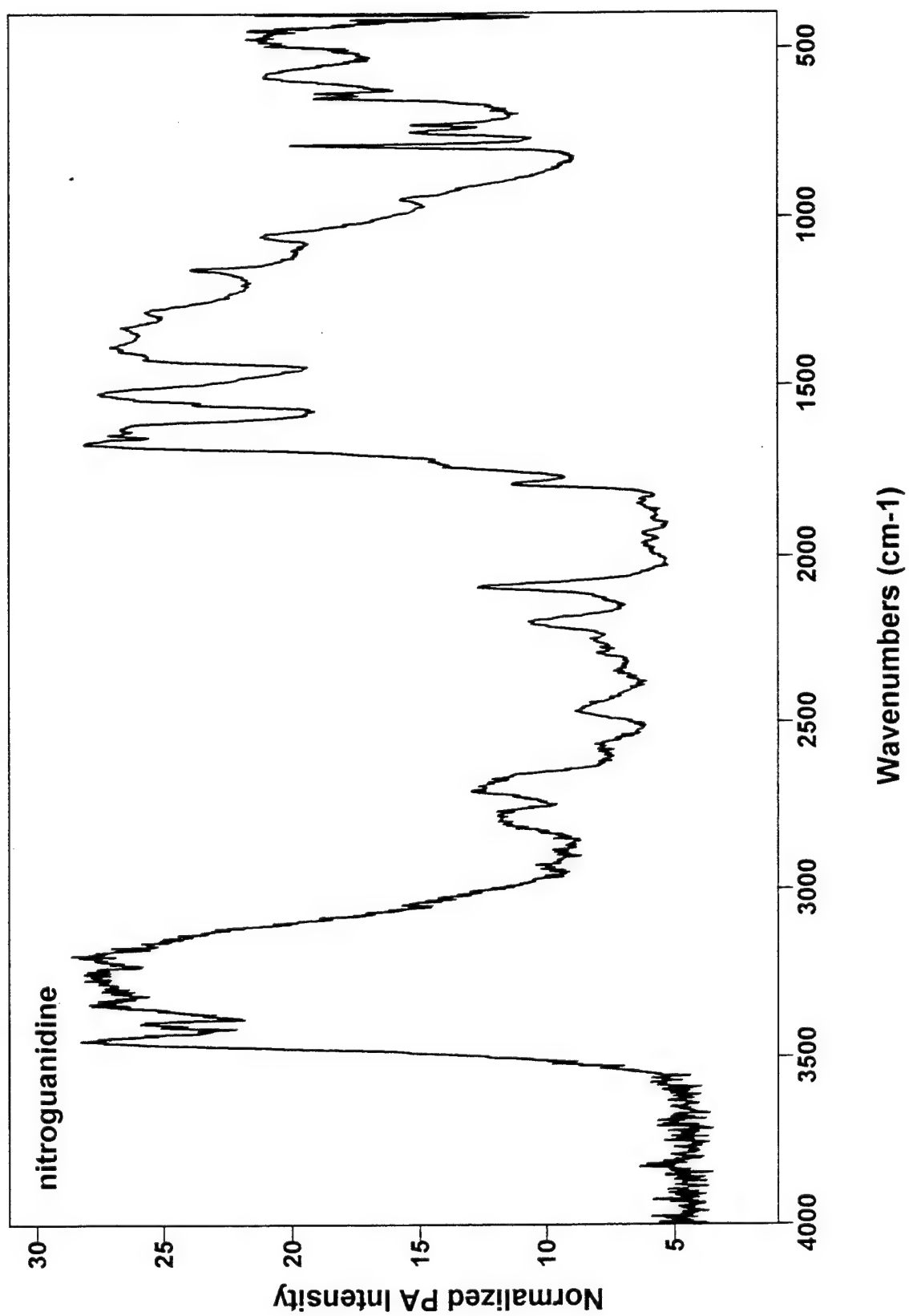



APPENDIX C:
SPECTRA OF BINDERS


INTENTIONALLY LEFT BLANK.

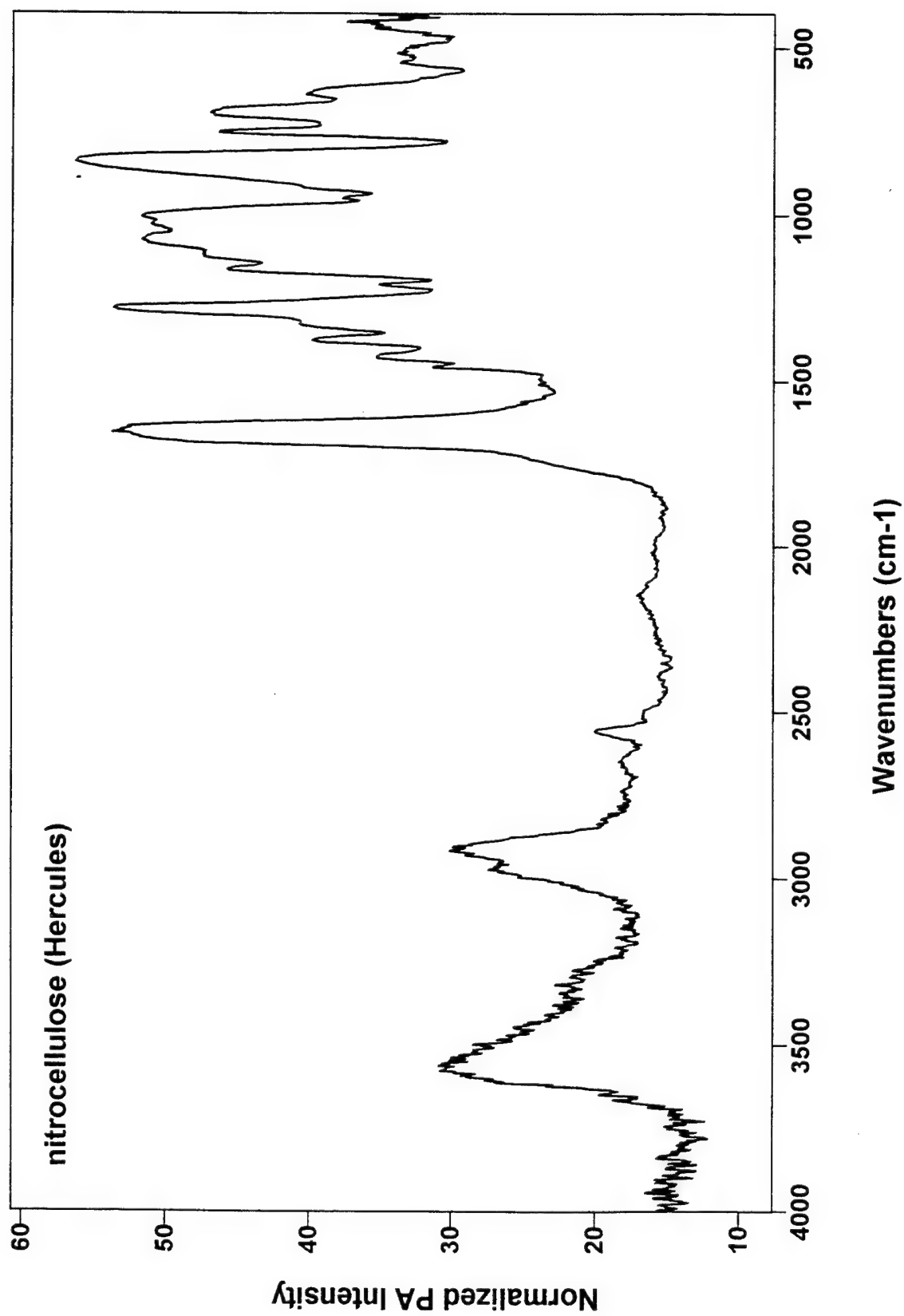


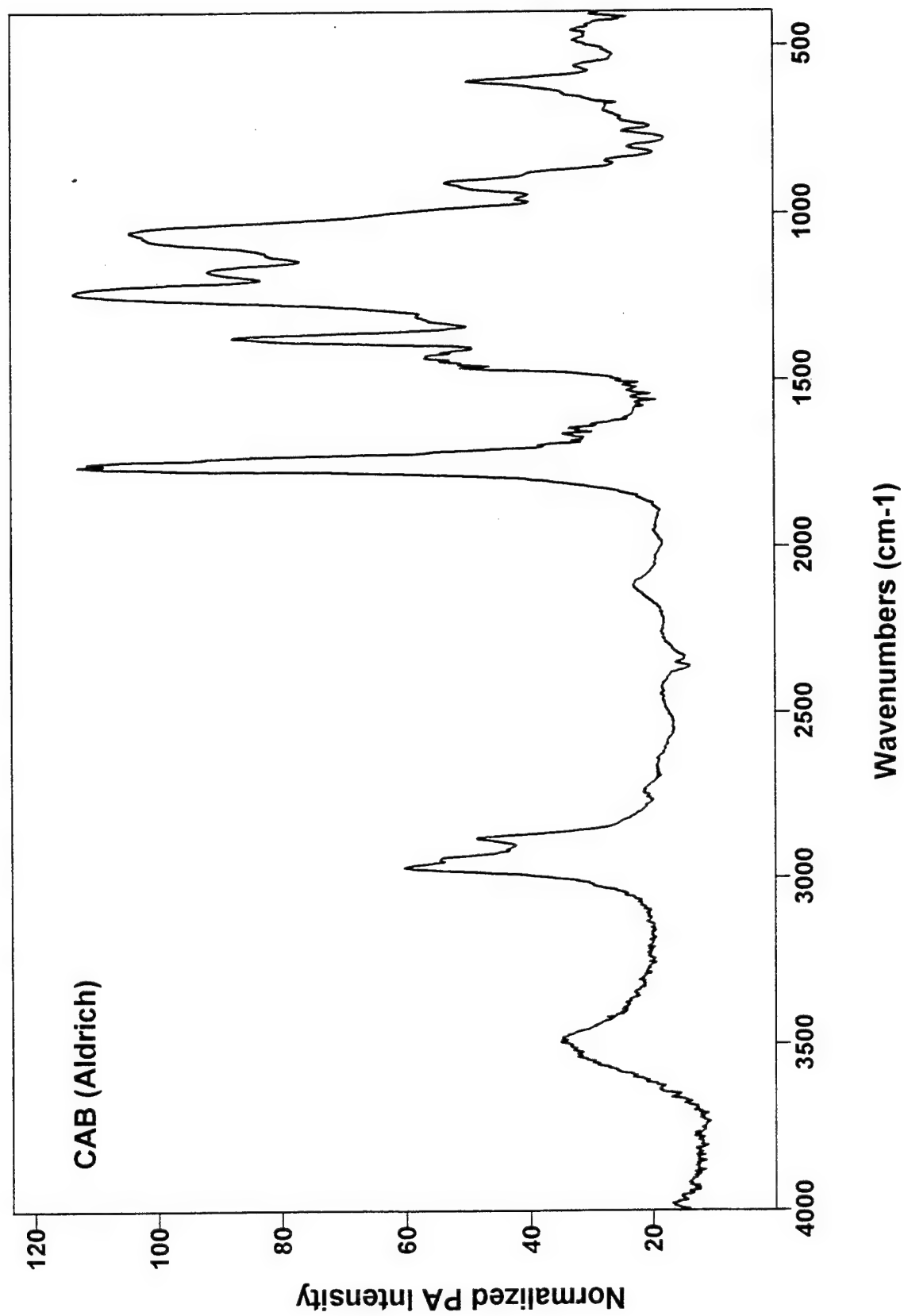


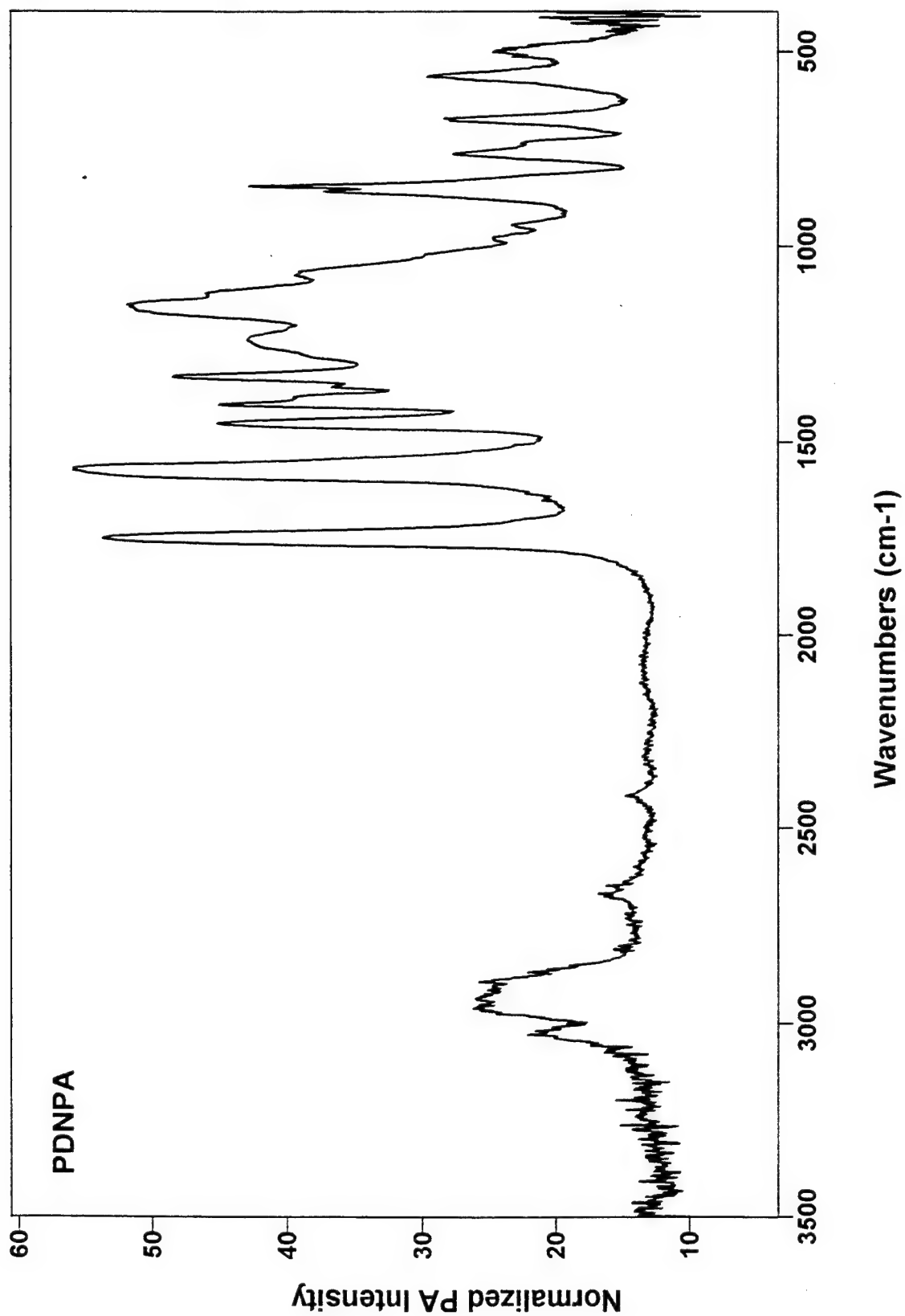


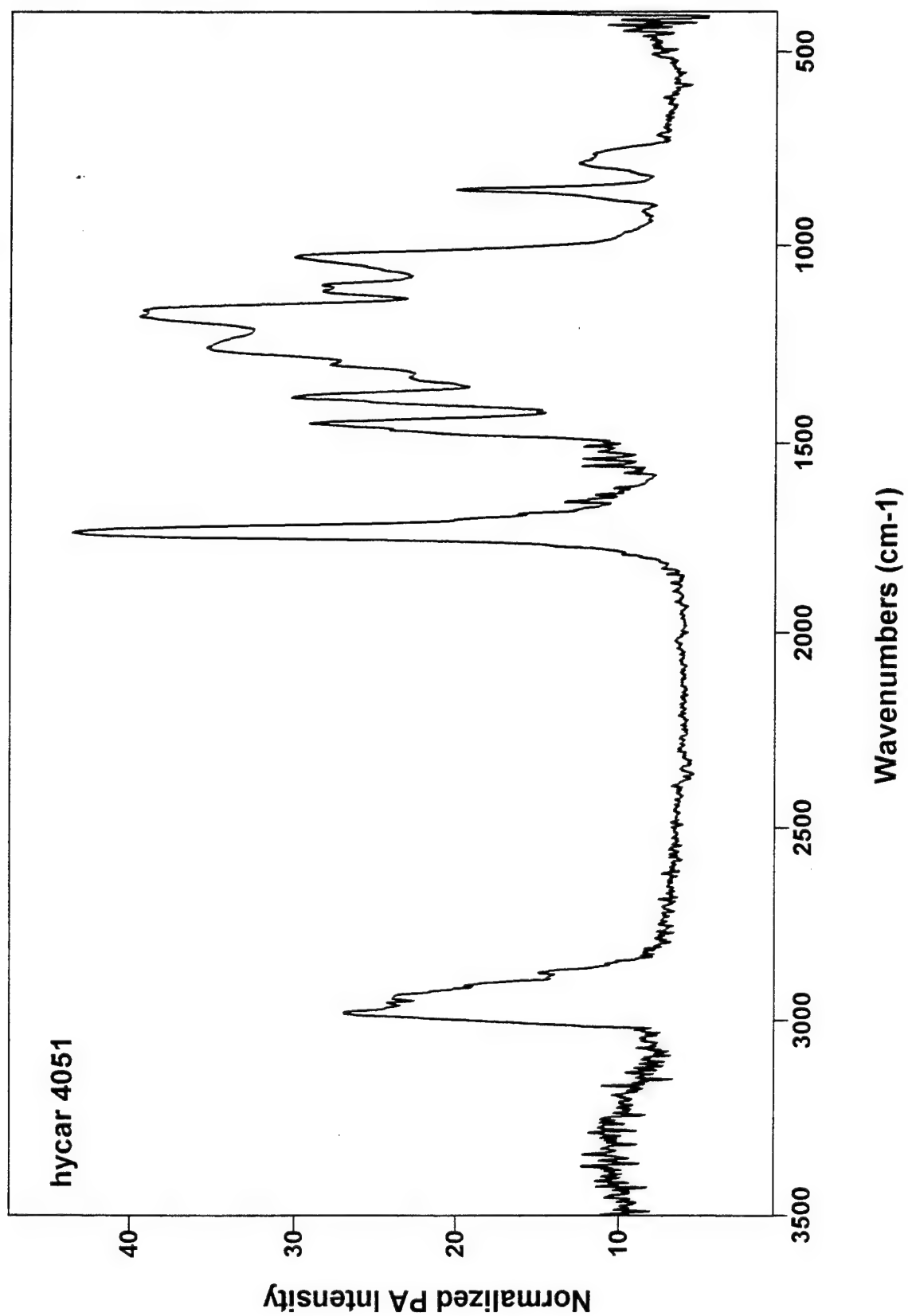


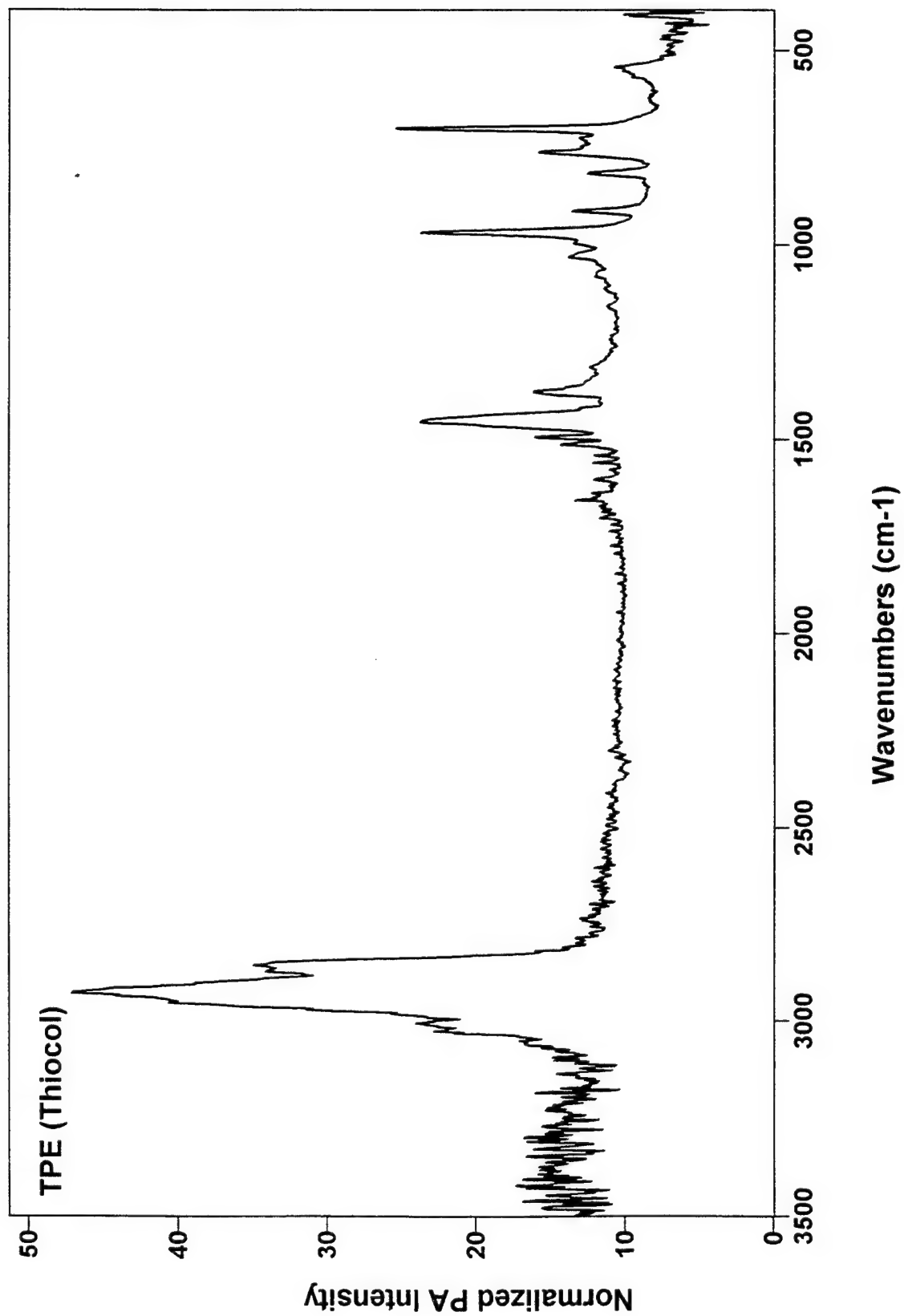

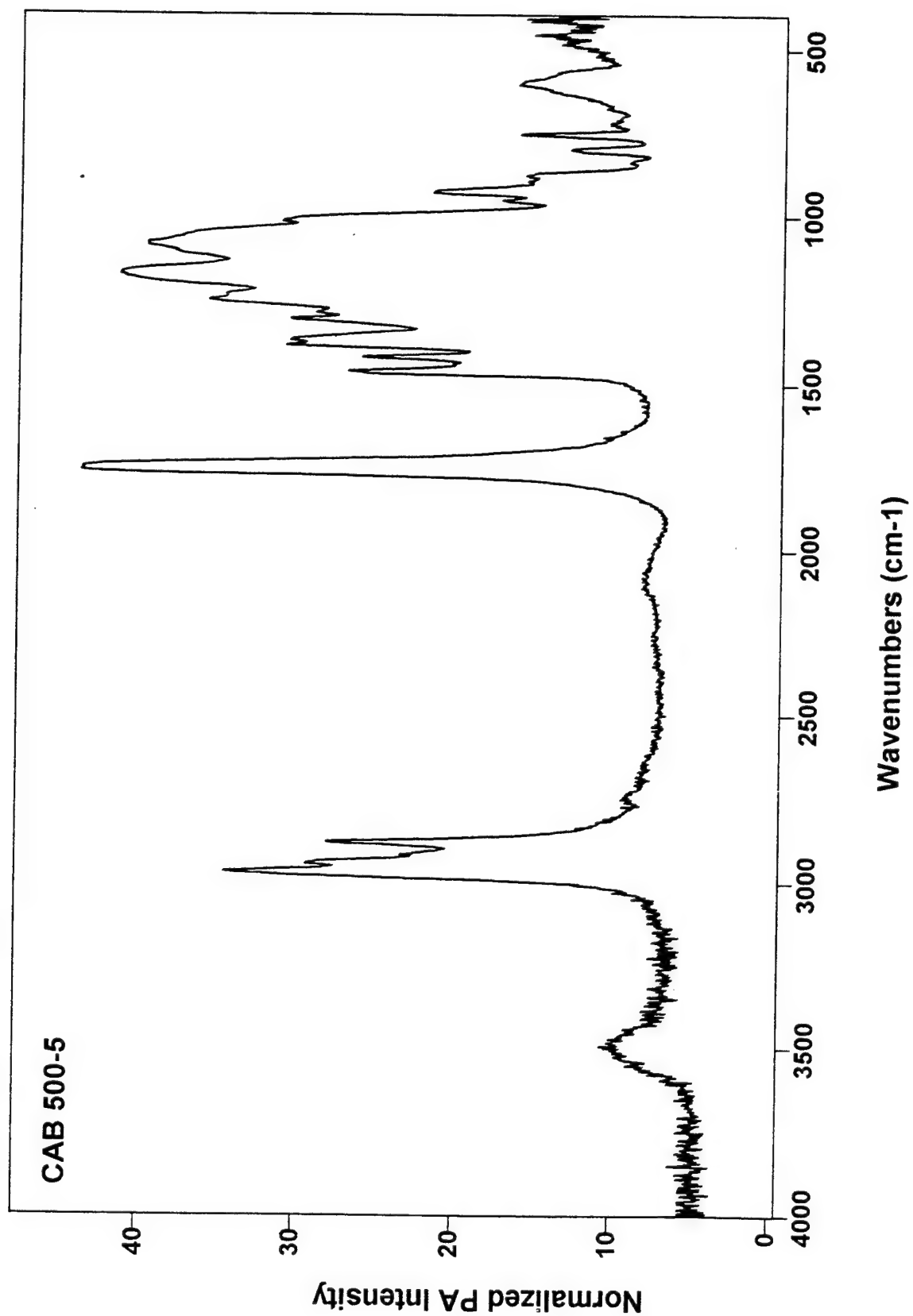

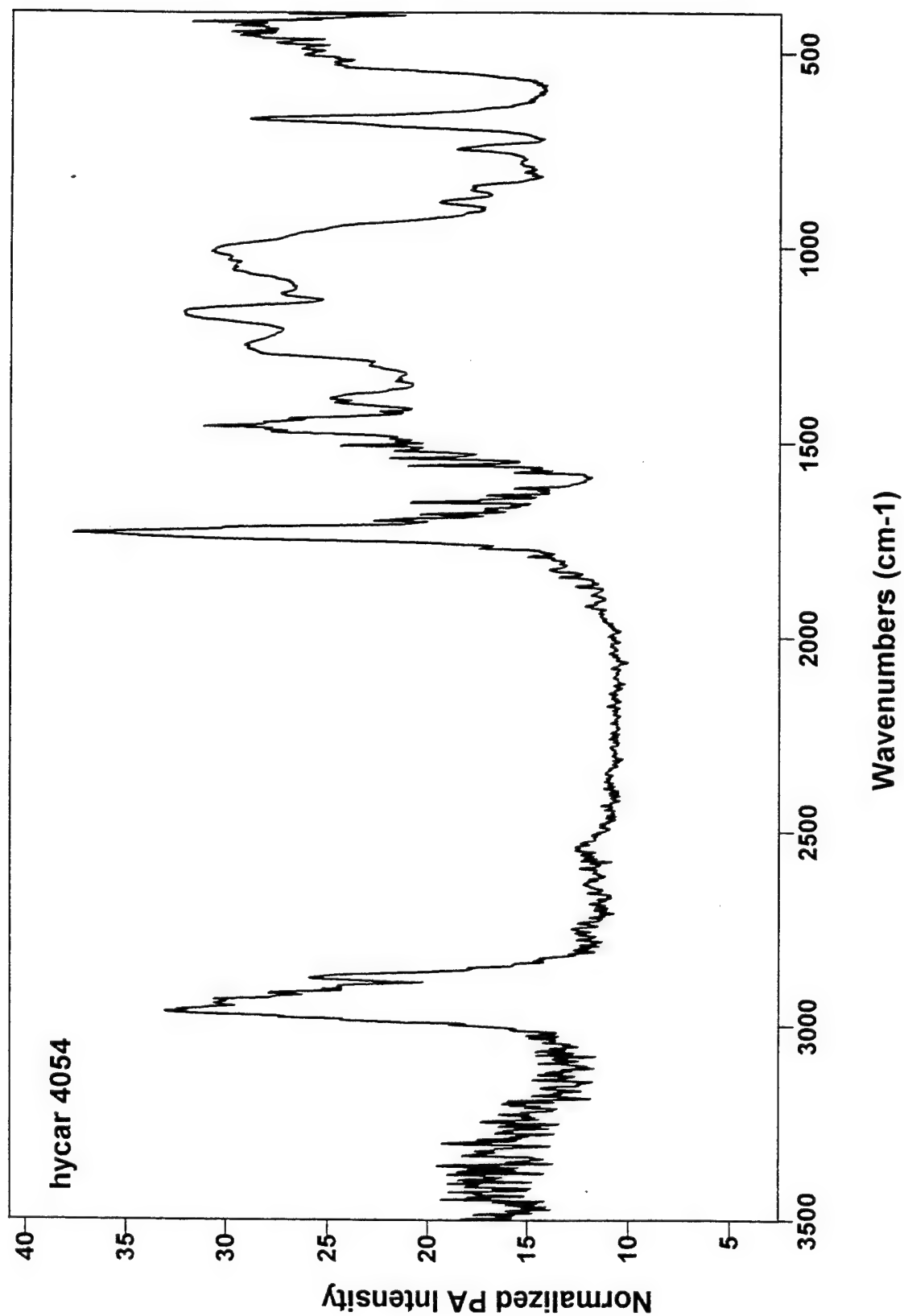

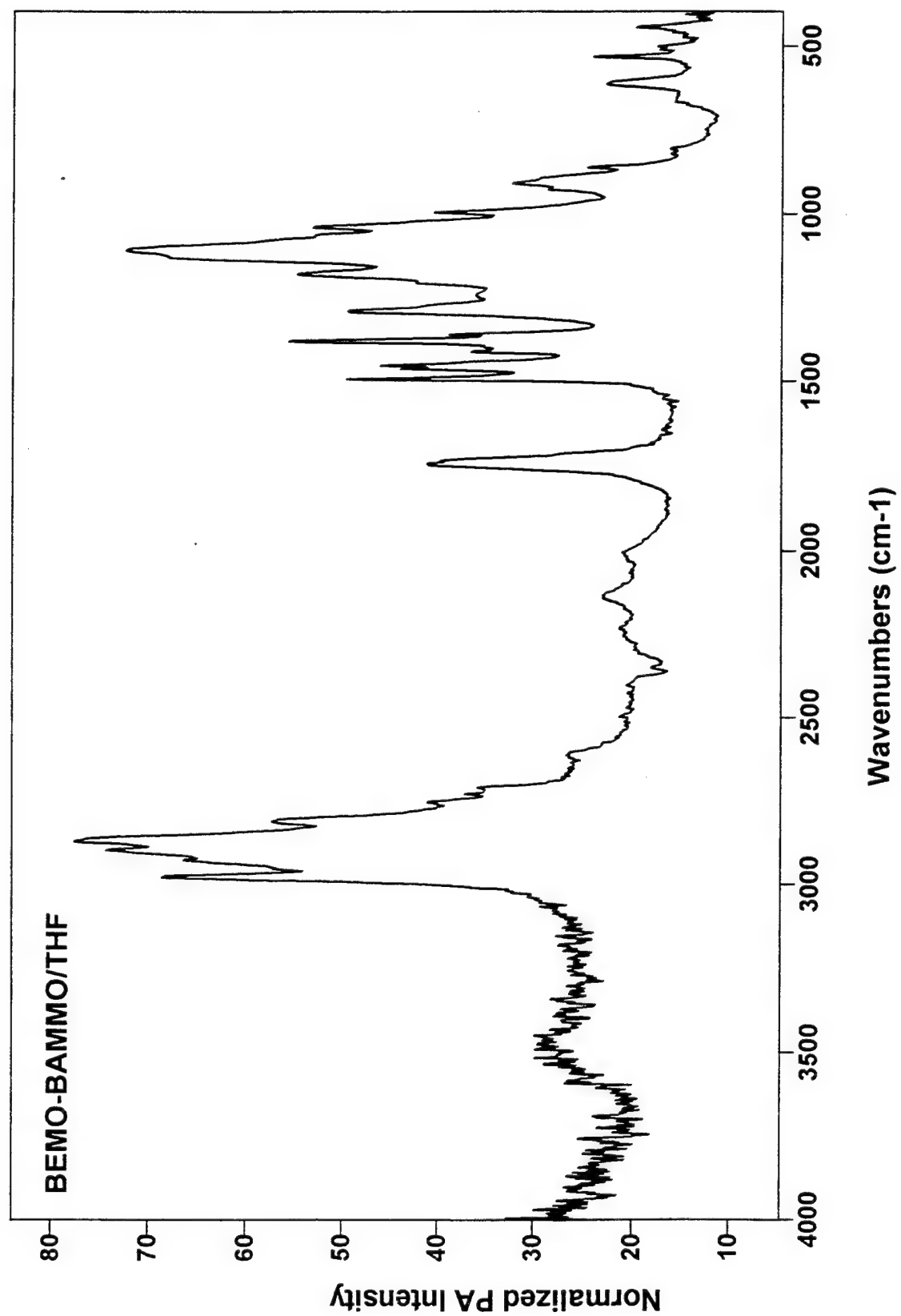


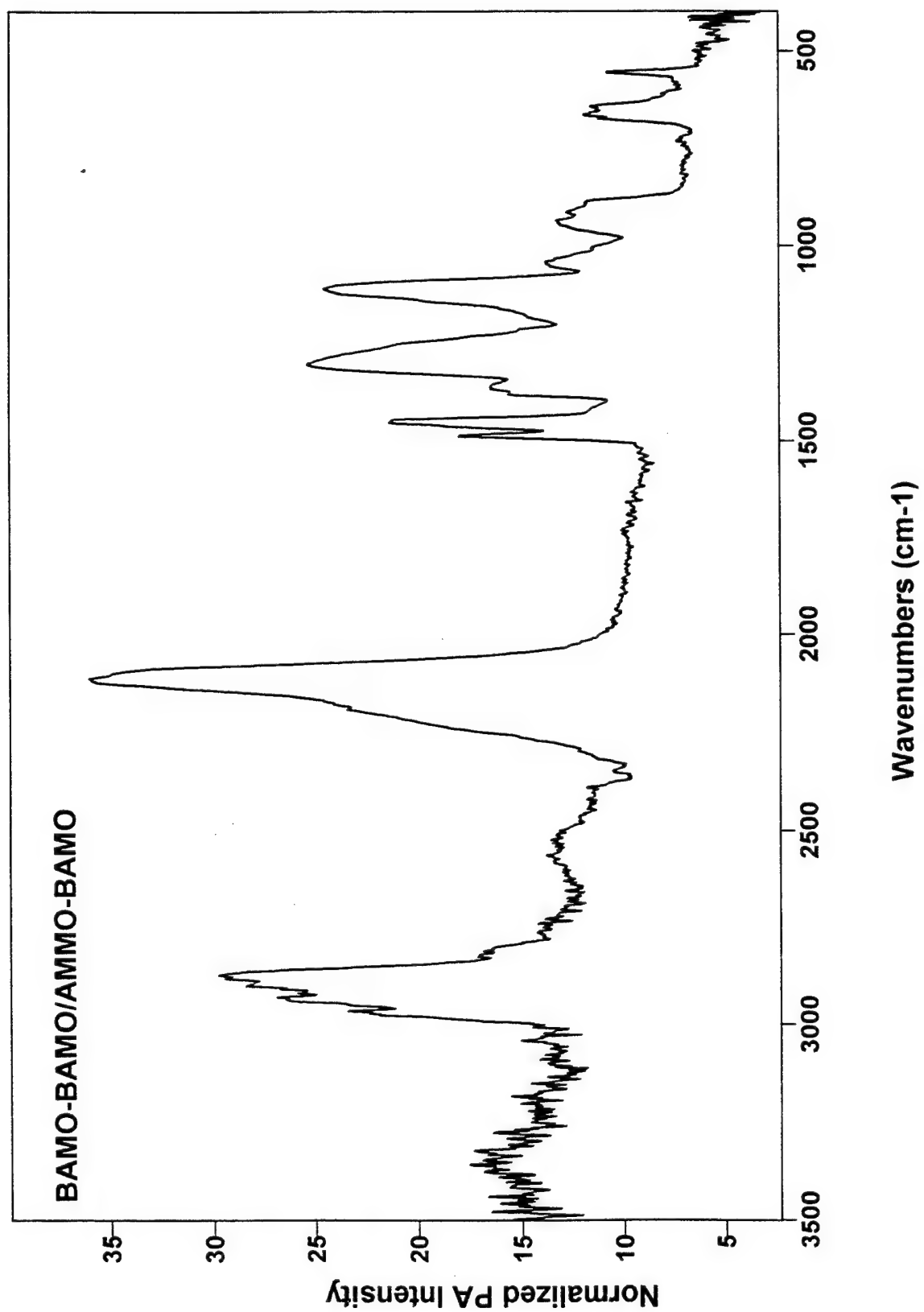



INTENTIONALLY LEFT BLANK.


APPENDIX D:
SPECTRA OF PLASTICIZERS


INTENTIONALLY LEFT BLANK.

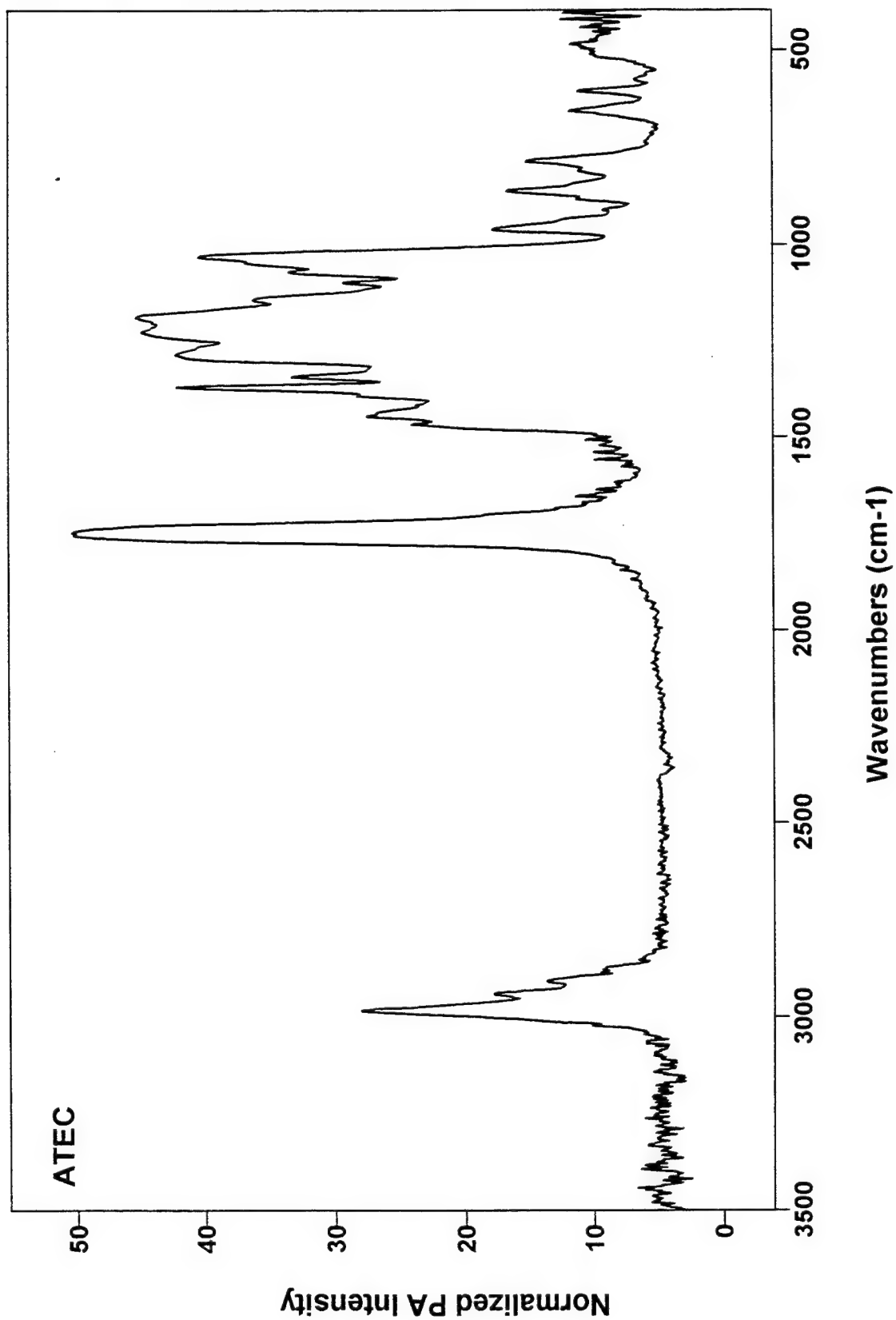

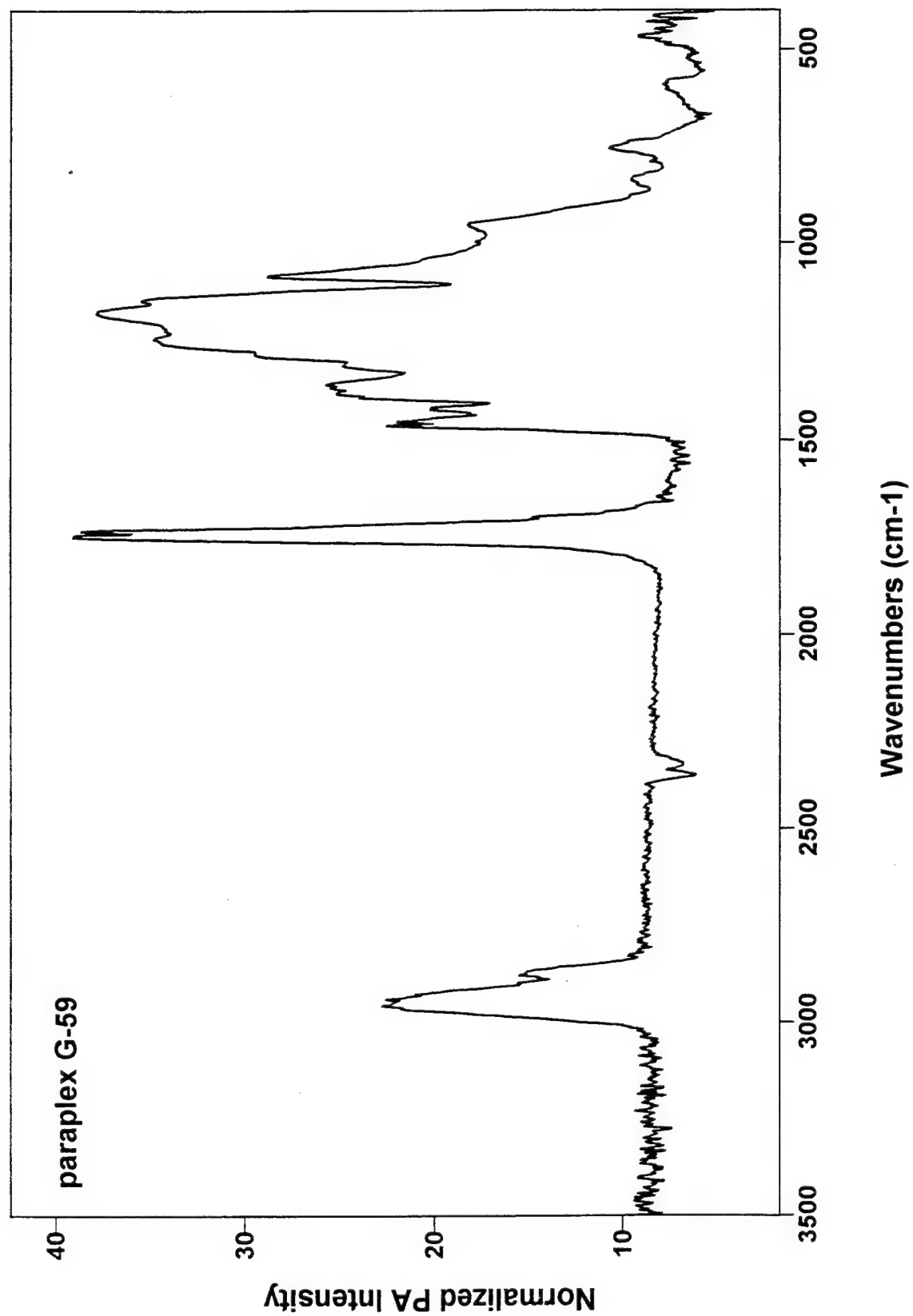

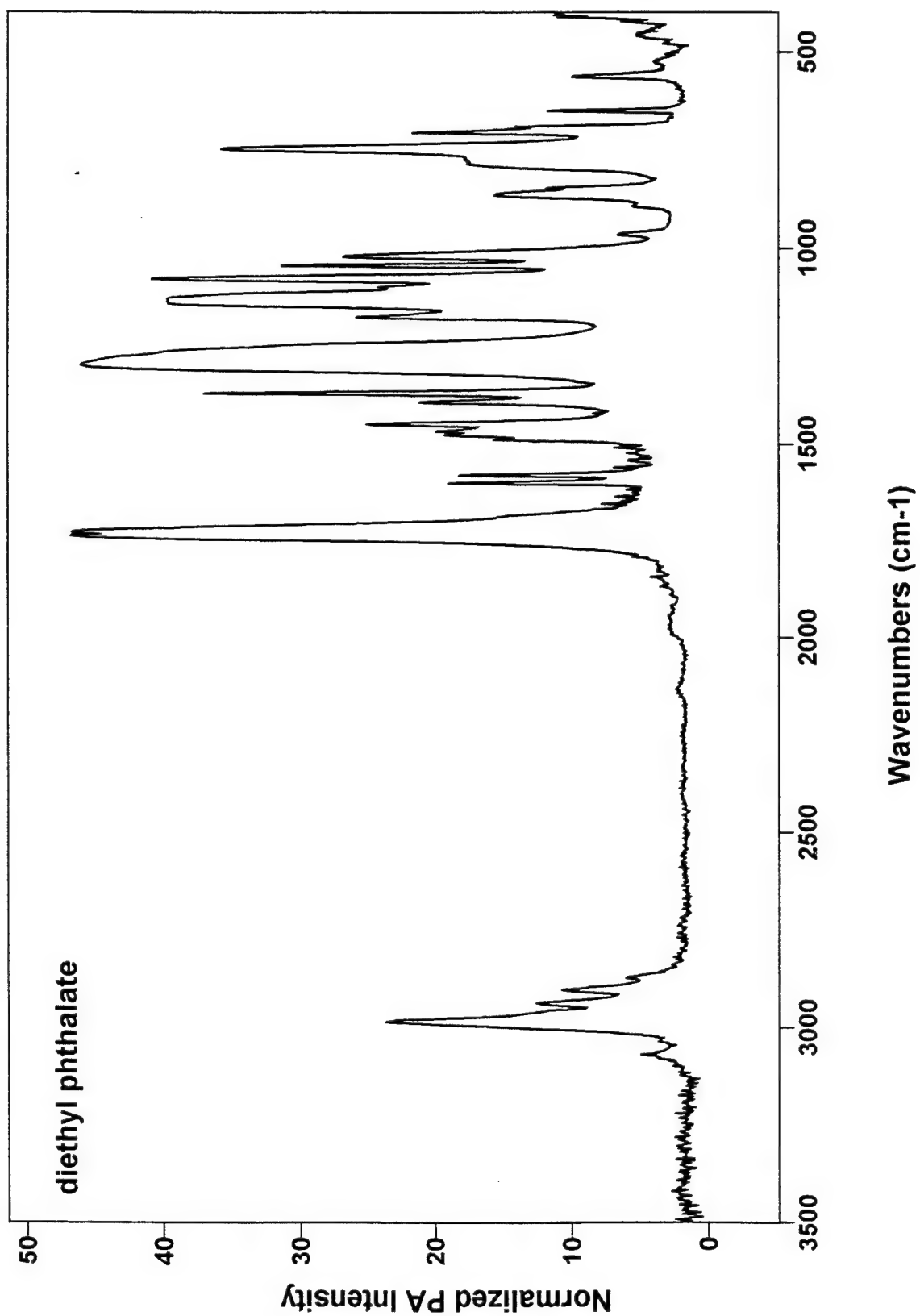

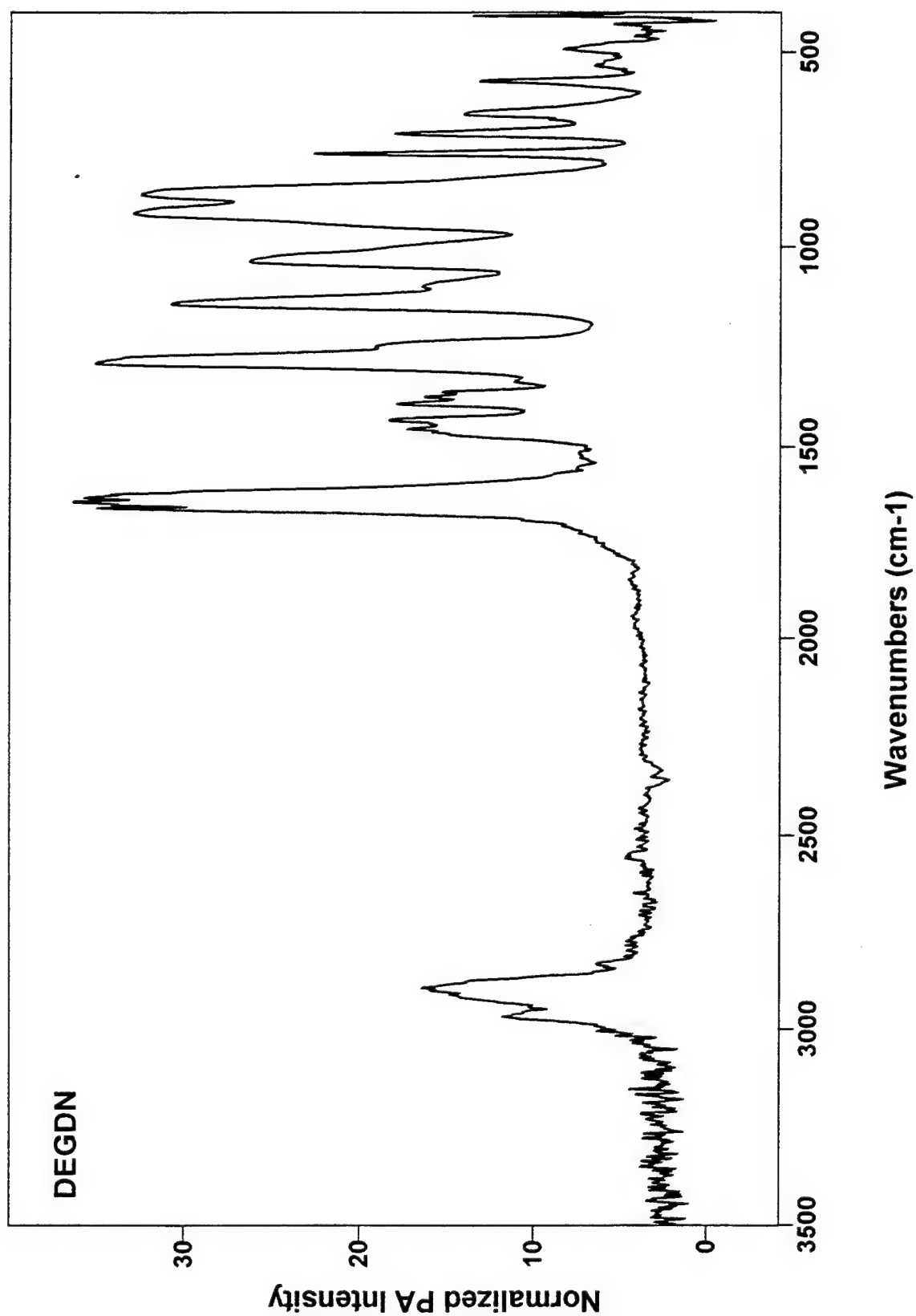


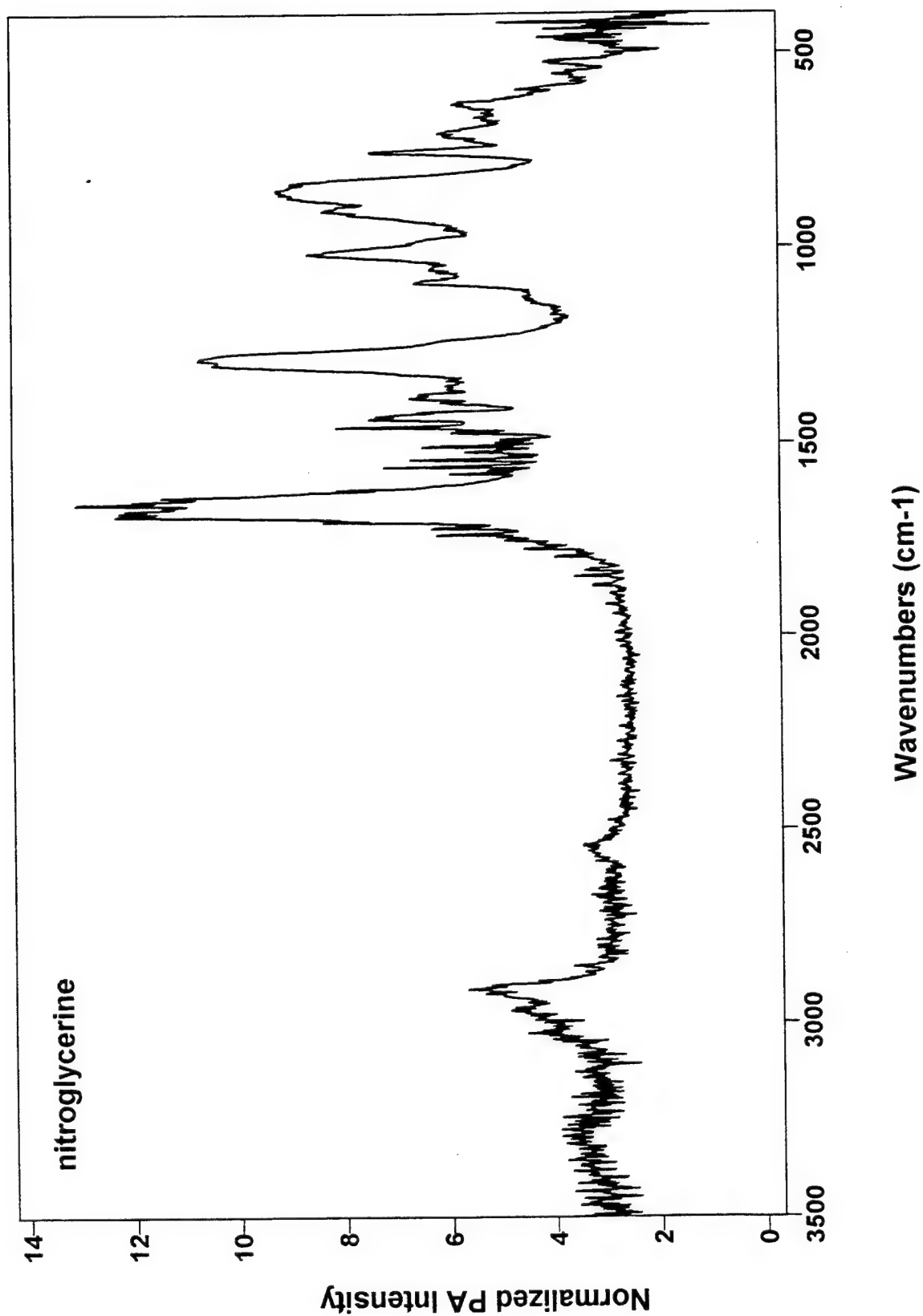



INTENTIONALLY LEFT BLANK.

APPENDIX E:
SPECTRA OF STABILIZERS

INTENTIONALLY LEFT BLANK.

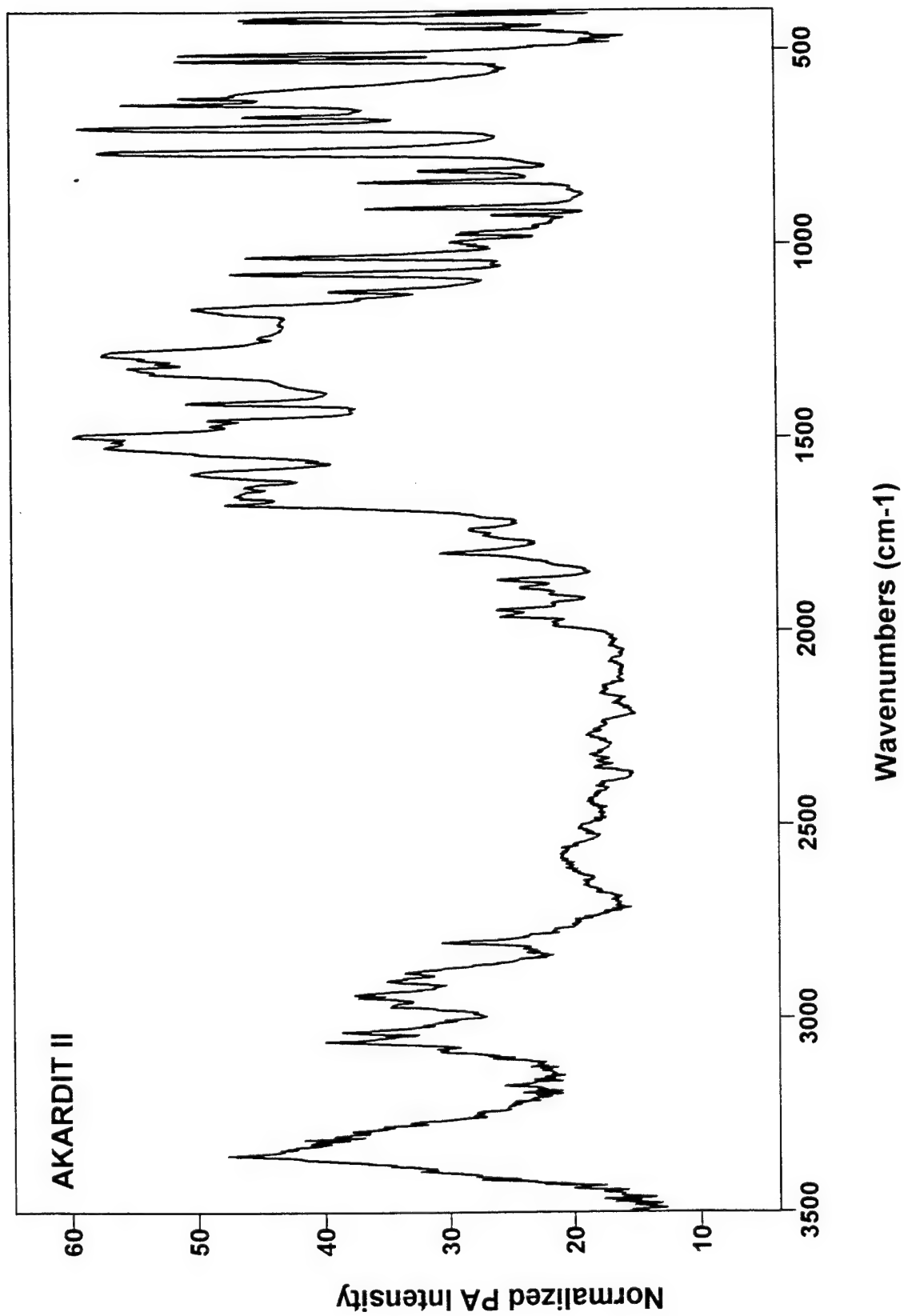

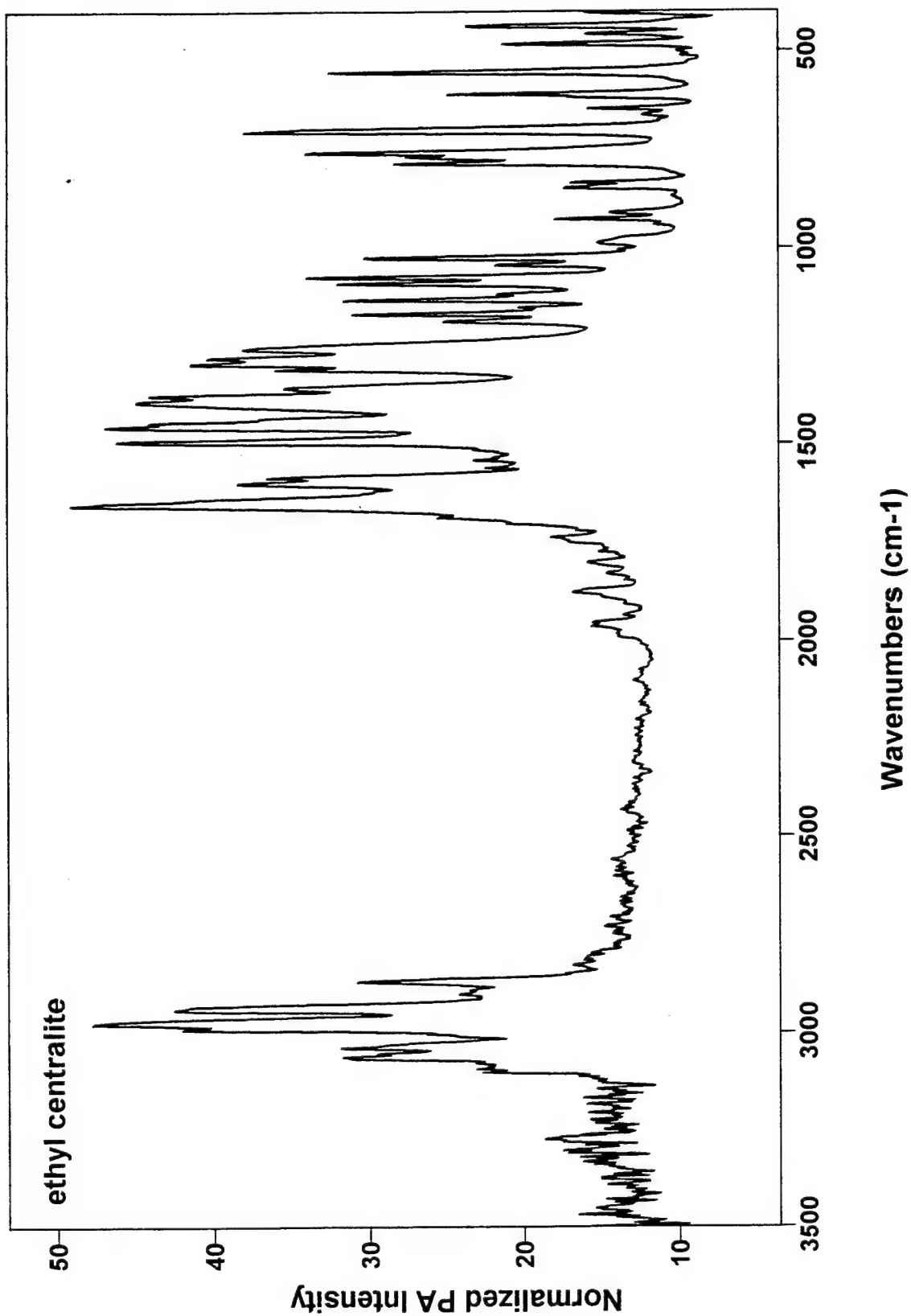


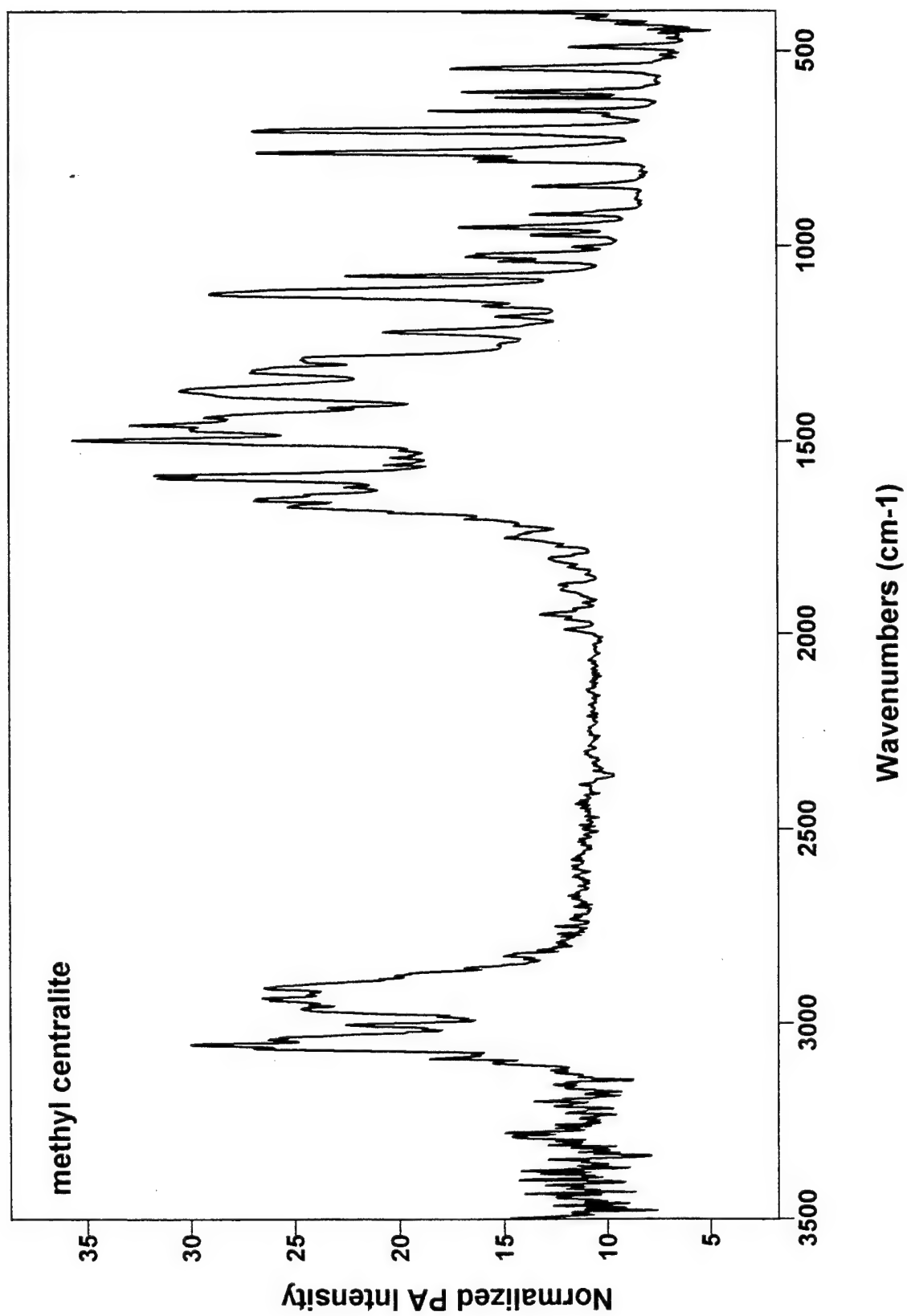



INTENTIONALLY LEFT BLANK.


APPENDIX F:
SPECTRA OF OTHERS


INTENTIONALLY LEFT BLANK.

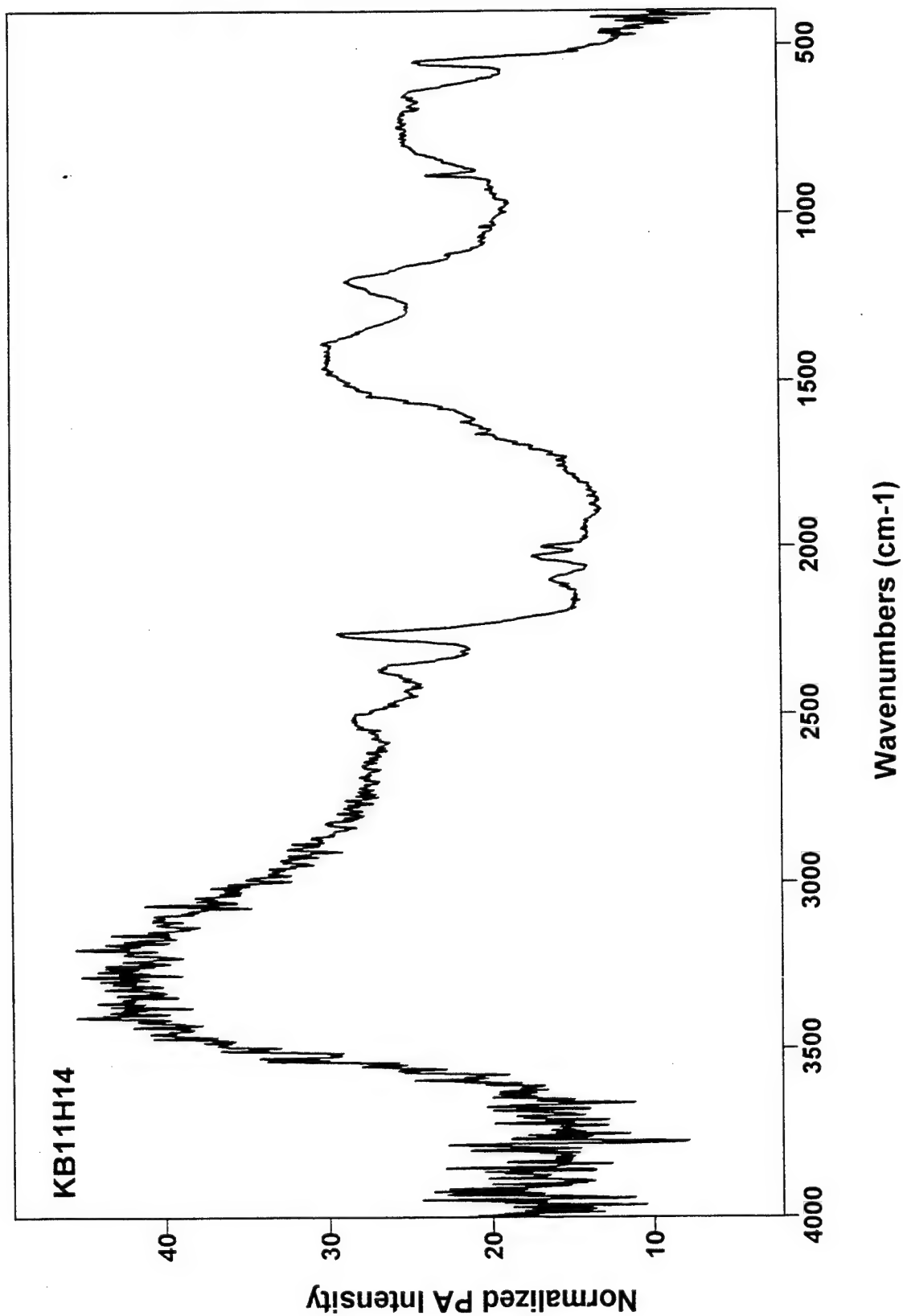


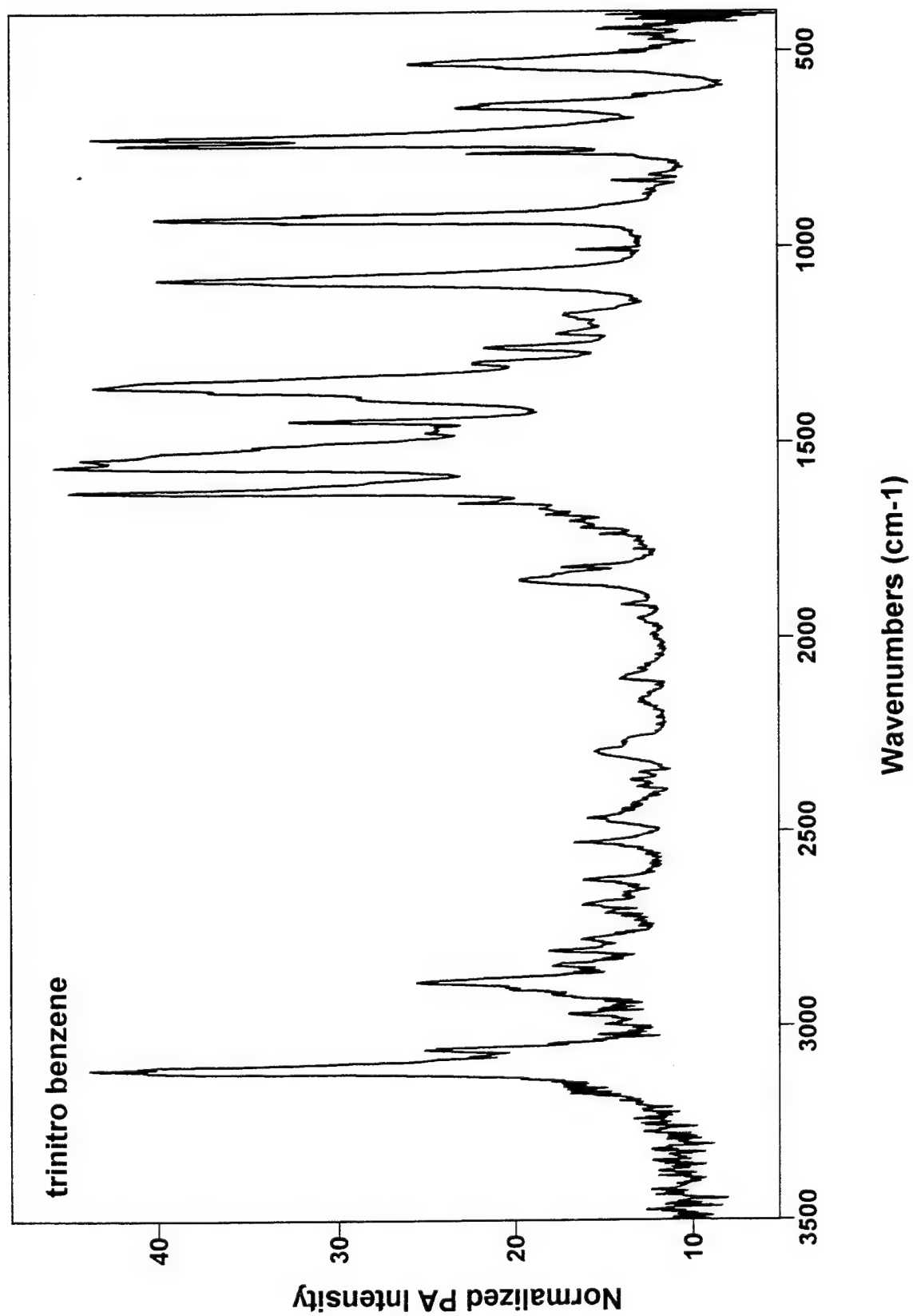


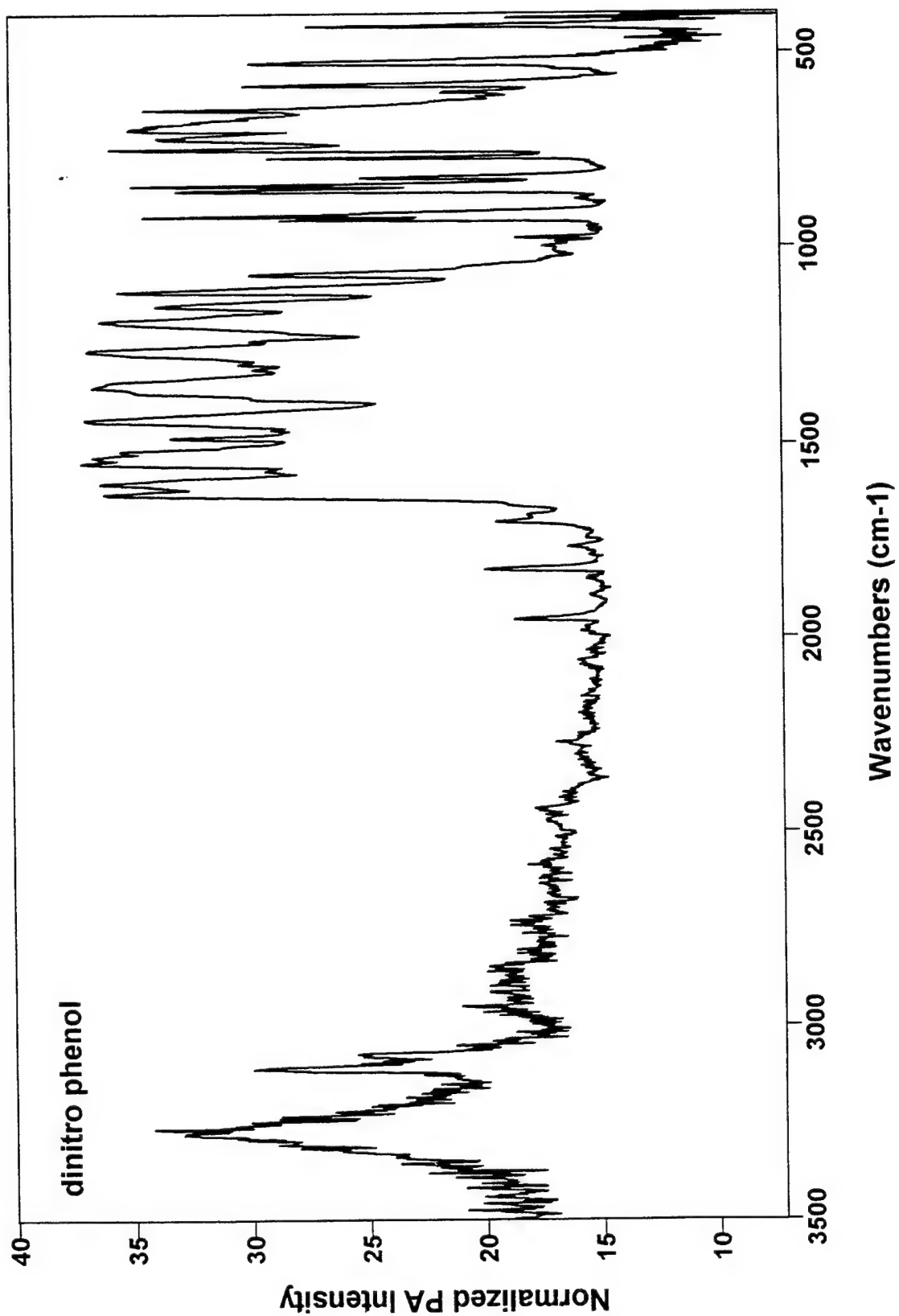

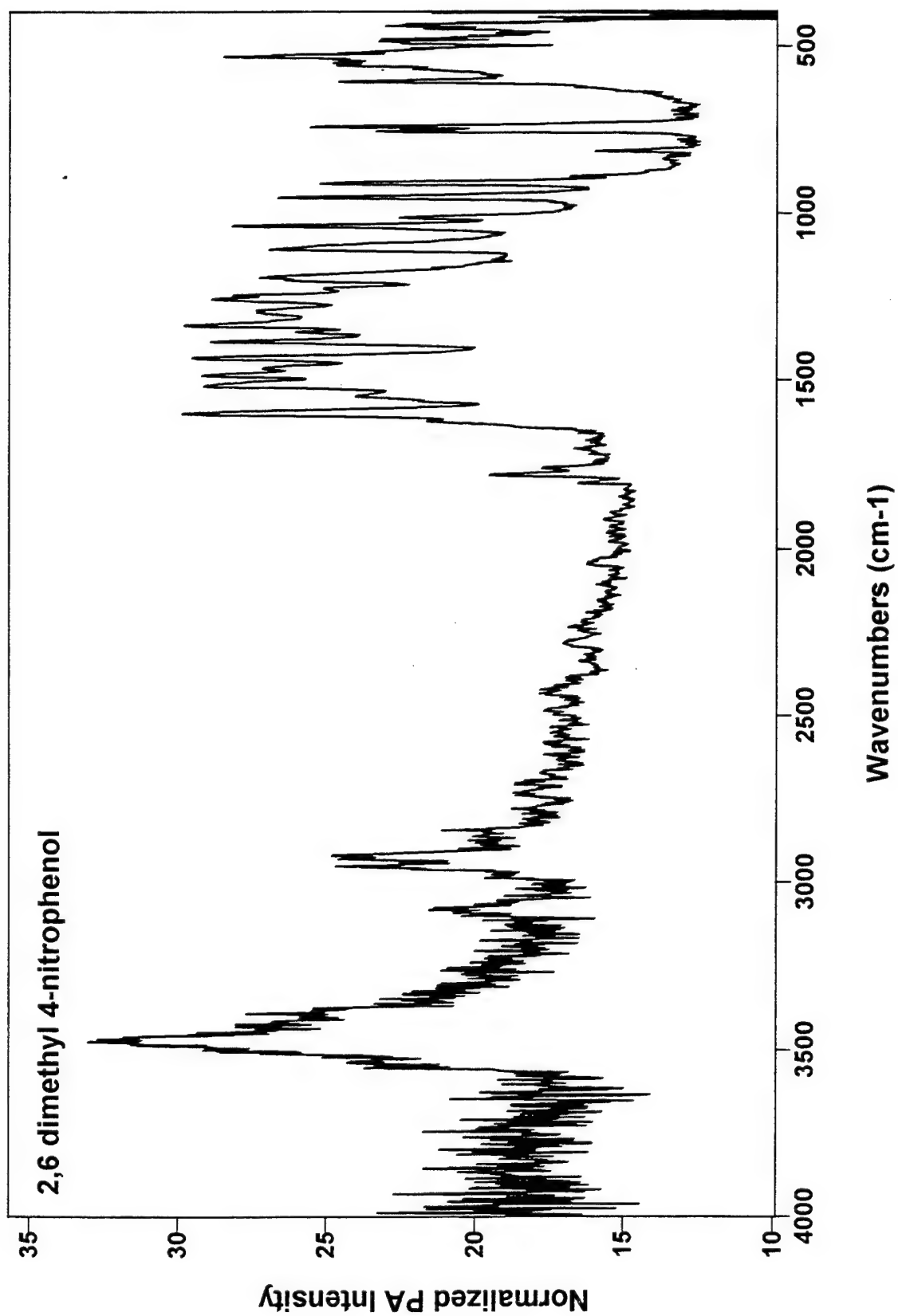

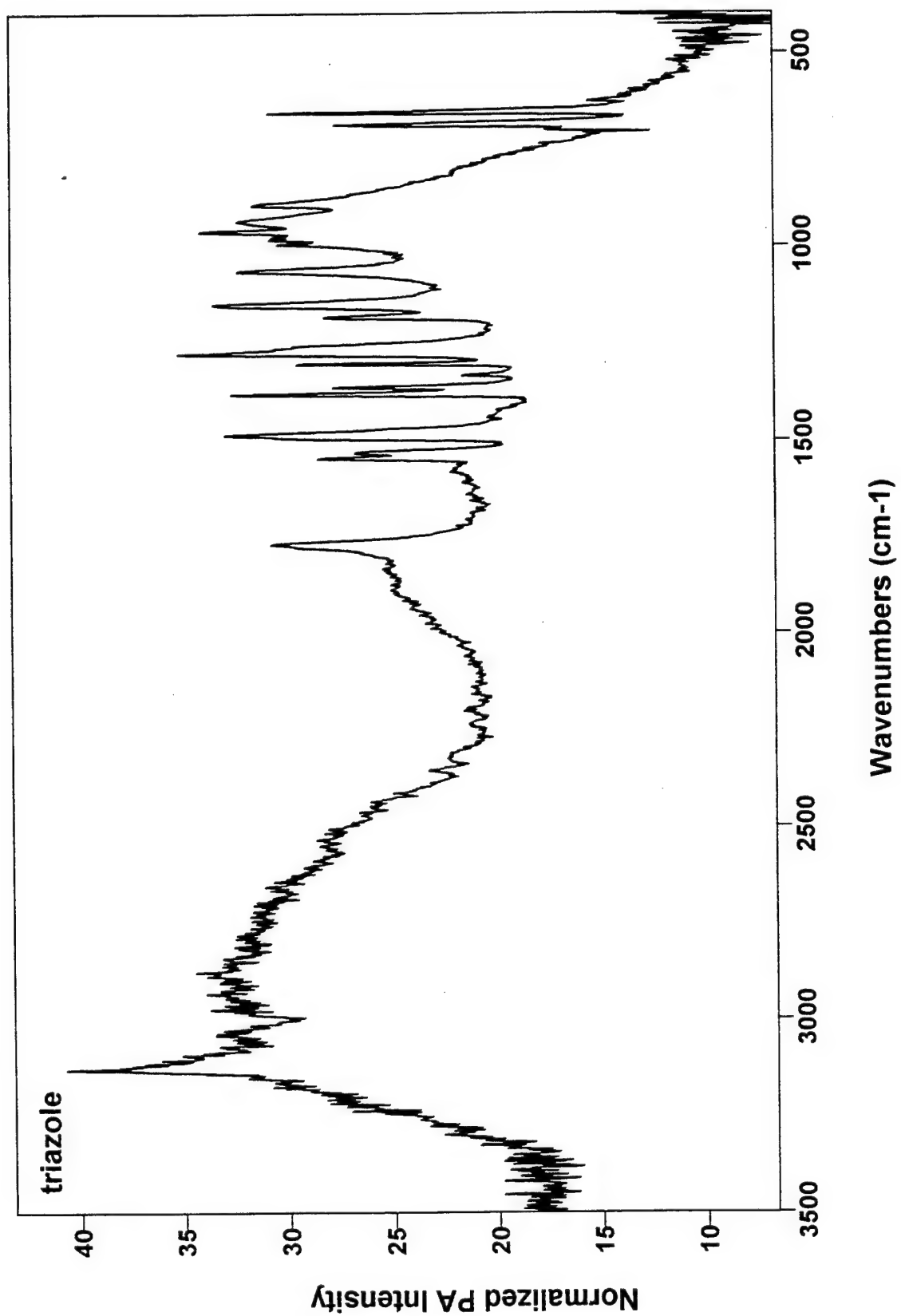

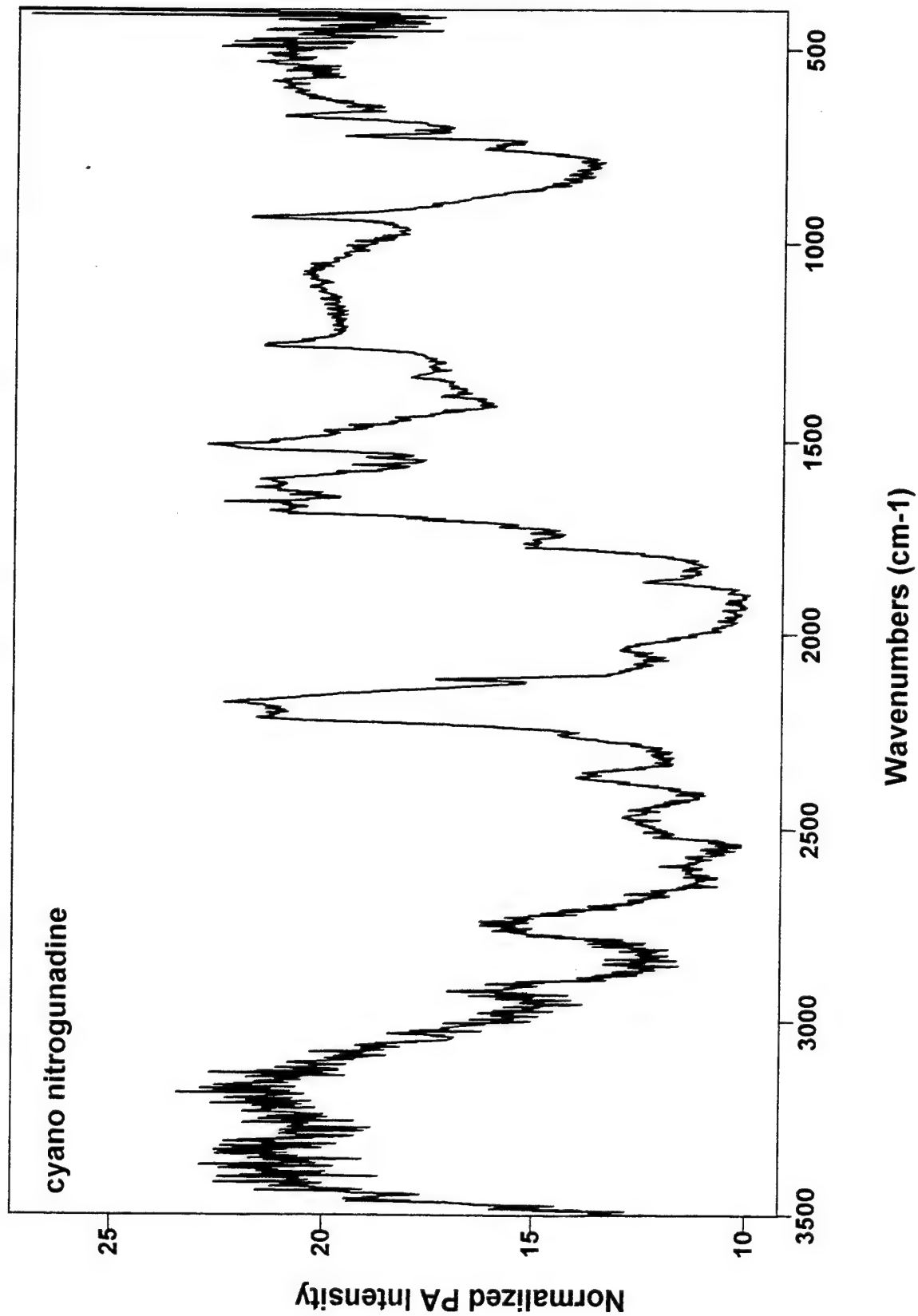

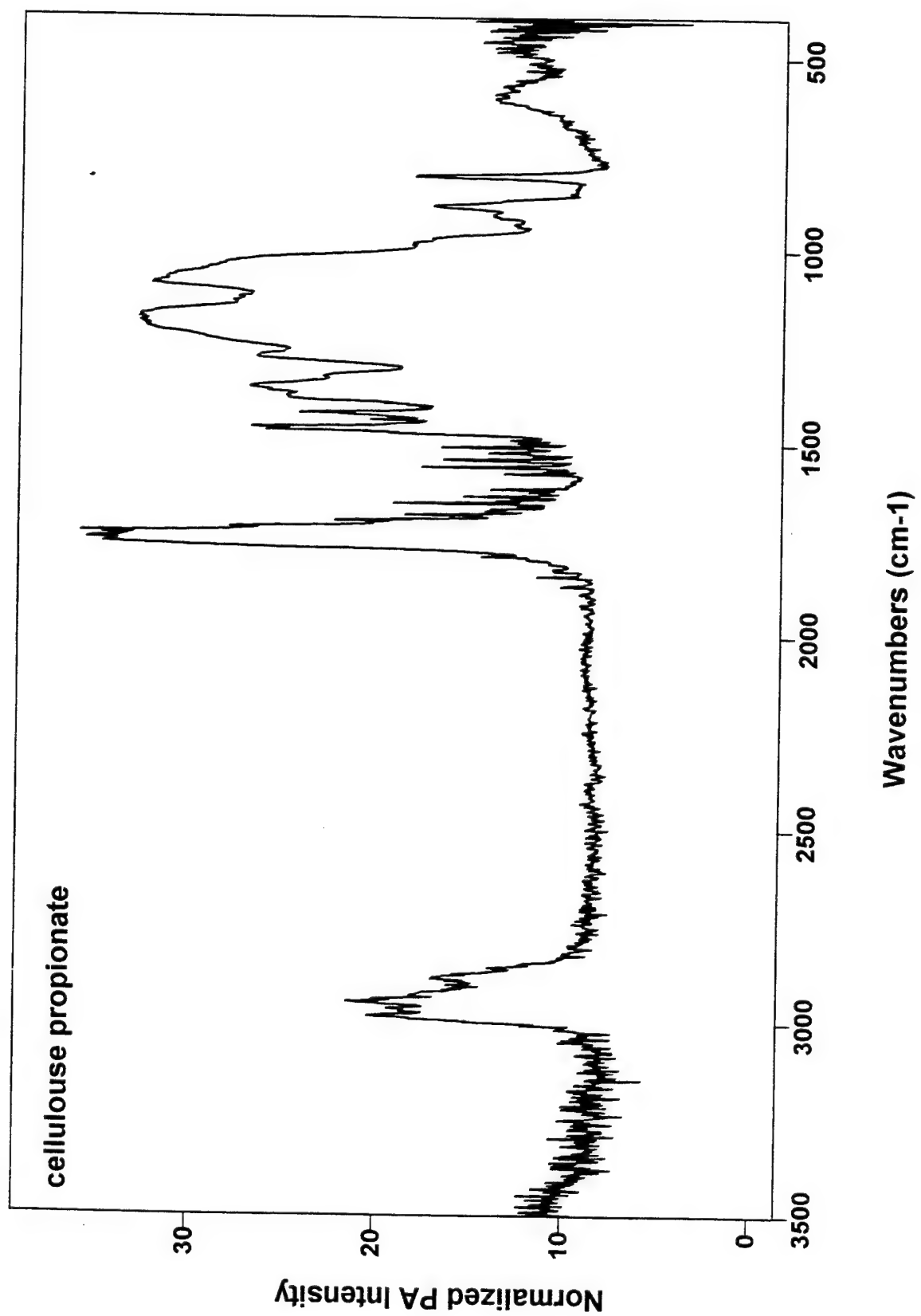

INTENTIONALLY LEFT BLANK.

LIST OF ABBREVIATIONS

AKARDIT II	-	Methyl Diphenyl Urea
AMMO	-	Azido Methyl Methyl Oxitane
ATEC	-	Acetyl Triethyl Citrate
ATR	-	Attenuated Total Reflectance
BAMO	-	Bis Azido Methyl Oxitane
BDNPA/F	-	Bis-Dinitro Propyl Acetal/Formal
BEMO	-	Bis Ethyl Methyl Oxitane
CAB	-	Cellulose Acetate Butyrate
CNNQ	-	Cyanonitroguanadine
CP	-	Cellulose Propionate
DANPE	-	Diazodo-Nitrazo-Pentane
DEGDN	-	Diethylene Glycol Dinitrate
DMNP	-	2,6, Dimethyl 4-Nitrophenol
DNP	-	Dinitrophenol
DPA	-	Diphenyl Amine
EC	-	Ethyl Centralite
FTIR	-	Fourier Transform Infrared
HELP	-	High Energy Low Vulnerability Propellant
HMX	-	High Melting Explosive
IR	-	Infrared
JAG	-	JA2 with Paraplex G-59
JAX	-	JA2 with RDX
KS	-	Potassium Sulfate
MC	-	Methyl Centralite
NC	-	Nitrocellulose
NG	-	Nitroglycerine
NQ	-	Nitroguanidine
PA	-	Photoacoustic
PDNPA	-	Poly (Dinitro Propyl Acrylate)
RDX	-	Royal Developmental Explosive
TNAZ	-	1,3,3-Trinitroazetidine
TNB	-	Trinitobenzene
TNT	-	Trinitrotoluene
TPE	-	Thermoplastic Elastomer

INTENTIONALLY LEFT BLANK.

<u>NO. OF COPIES</u>	<u>ORGANIZATION</u>
2	DEFENSE TECHNICAL INFO CTR ATTN DTIC DDA 8725 JOHN J KINGMAN RD STE 0944 FT BELVOIR VA 22060-6218
1	DIRECTOR US ARMY RESEARCH LAB ATTN AMSRL CS AL TP 2800 POWDER MILL RD ADELPHI MD 20783-1145
1	DIRECTOR US ARMY RESEARCH LAB ATTN AMSRL CS AL TA 2800 POWDER MILL RD ADELPHI MD 20783-1145
3	DIRECTOR US ARMY RESEARCH LAB ATTN AMSRL CI LL 2800 POWDER MILL RD ADELPHI MD 20783-1145
	<u>ABERDEEN PROVING GROUND</u>
2	DIR USARL ATTN AMSRL CI LP (305)

<u>NO. OF COPIES</u>	<u>ORGANIZATION</u>
1	HQDA PENTAGON ATTN SARD TT DR F MILTON WASHINGTON DC 20310-0103
1	HQDA PENTAGON ATTN SARD TT MR J APPEL WASHINGTON DC 20310-0103
1	HQDA OASA RDA PENTAGON ROOM 3E486 ATTN DR C H CHURCH WASHINGTON DC 20310-0103
4	COMMANDER US ARMY RESEARCH OFFICE ATTN R GHIRARDELLI D MANN R SINGLETON R SHAW PO BOX 12211 RESEARCH TRIANGLE PARK NC 27709-2211
1	DIRECTOR ARMY RESEARCH OFFICE ATTN AMXRO RT IP LIB SERVICES PO BOX 12211 RESEARCH TRIANGLE PARK NC 27709-2211
2	COMMANDER US ARMY ARDEC ATTN SMCAR AEE B D S DOWNS PICATINNY ARSENAL NJ 07806-5000
2	COMMANDER US ARMY ARDEC ATTN SMCAR AEE J A LANNON PICATINNY ARSENAL NJ 07806-5000
1	COMMANDER US ARMY ARDEC ATTN SMCAR AEE BR L HARRIS PICATINNY ARSENAL NJ 07806-5000
2	COMMANDER US ARMY MISSILE COMMAND ATTN AMSMI RD PR E A R MAYKUT AMSMI RD PR P R BETTS REDSTONE ARSENAL AL 35809

<u>NO. OF COPIES</u>	<u>ORGANIZATION</u>
1	OFFICE OF NAVAL RESEARCH DEPARTMENT OF THE NAVY ATTN R S MILLER CODE 432 800 N QUINCY STREET ARLINGTON VA 22217
1	COMMANDER NAVAL AIR SYSTEMS COMMAND ATTN J RAMNARACE AIR 54111C WASHINGTON DC 20360
2	COMMANDER NAVAL SURFACE WARFARE CENTER ATTN R BERNECKER R 13 G B WILMOT R 16 SILVER SPRING MD 20903-5000
5	COMMANDER NAVAL RESEARCH LABORATORY ATTN M C LIN J MCDONALD E ORAN J SHNUR R J DOYLE CODE 6110 WASHINGTON DC 20375
2	COMMANDER NAVAL WEAPONS CENTER ATTN T BOGGS CODE 388 T PARR CODE 3895 CHINA LAKE CA 93555-6001
1	SUPERINTENDENT NAVAL POSTGRADUATE SCHOOL DEPT OF AERONAUTICS ATTN D W NETZER MONTEREY CA 93940
3	AL LSCF ATTN R CORLEY R GEISLER J LEVINE EDWARDS AFB CA 93523-5000
1	AFOSR BOLLING AIR FORCE BASE ATTN J M TISHKOFF WASHINGTON DC 20332

NO. OF
COPIES ORGANIZATION

1 OSD SDIO IST
PENTAGON
ATTN L CAVENY
WASHINGTON DC 20301-7100

1 COMMANDANT
USAFAS
ATTN ATSF TSM CN
FORT SILL OK 73503-5600

1 UNIV OF DAYTON RSCH INSTITUTE
ATTN D CAMPBELL
AL PAP
EDWARDS AFB CA 93523

1 NASA
LANGLEY RESEARCH CENTER
ATTN G B NORTHAM MS 168
LANGLEY STATION
HAMPTON VA 23365

4 NATIONAL BUREAU OF STANDARDS
US DEPARTMENT OF COMMERCE
ATTN J HASTIE
M JACOX
T KASHIWAGI
H SEMERJIAN
WASHINGTON DC 20234

2 DIRECTOR
LAWRENCE LIVERMORE NATIONAL LAB
ATTN C WESTBROOK
W TAO MS L 282
P O BOX 808
LIVERMORE CA 94550

1 DIRECTOR
LOS ALAMOS NATIONAL LAB
ATTN B NICHOLS T7 MS B284
P O BOX 1663
LOS ALAMOS NM 87545

2 PRINCETON COMBUSTION
RESEARCH LABORATORIES INC
ATTN N A MESSINA
M SUMMERFIELD
PRINCETON CORPORATE PLAZA
BLDG IV SUITE 119
11 DEERPARK DRIVE
MONMOUTH JUNCTION NJ 08852

NO. OF
COPIES ORGANIZATION

3 DIRECTOR
SANDIA NATIONAL LABORATORIES
DIVISION 8354
ATTN S JOHNSTON
P MATTERN
D STEPHENSON
LIVERMORE CA 94550

1 BRIGHAM YOUNG UNIVERSITY
DEPT OF CHEMICAL ENGINEERING
ATTN M W BECKSTEAD
PROVO UT 84058

1 CALIFORNIA INSTITUTE OF TECH
JET PROPULSION LABORATORY
ATTN L STRAND MS 125 224
4800 OAK GROVE DRIVE
PASADENA CA 91109

1 CALIFORNIA INSTITUTE OF TECHNOLOGY
ATTN F E C CULICK MC 301 46
204 KARMAN LAB
PASADENA CA 91125

1 UNIVERSITY OF CALIFORNIA
LOS ALAMOS SCIENTIFIC LAB
P O BOX 1663 MAIL STOP B216
LOS ALAMOS NM 87545

1 UNIVERSITY OF CALIFORNIA BERKELEY
CHEMISTRY DEPARMENT
ATTN C BRADLEY MOORE
211 LEWIS HALL
BERKELEY CA 94720

1 UNIVERSITY OF CALIFORNIA SAN DIEGO
ATTN F A WILLIAMS
AMES B010
LA JOLLA CA 92093

2 UNIV OF CALIFORNIA SANTA BARBARA
QUANTUM INSTITUTE
ATTN K SCHOFIELD
M STEINBERG
SANTA BARBARA CA 93106

1 UNIV OF COLORADO AT BOULDER
ENGINEERING CENTER
ATTN J DAILY
CAMPUS BOX 427
BOULDER CO 80309-0427

NO. OF COPIES	ORGANIZATION
3	UNIV OF SOUTHERN CALIFORNIA DEPT OF CHEMISTRY ATTN R BEAUDET S BENSON C WITTIG LOS ANGELES CA 90007
1	CORNELL UNIVERSITY DEPARTMENT OF CHEMISTRY ATTN T A COOL BAKER LABORATORY ITHACA NY 14853
1	UNIVERSITY OF DELAWARE CHEMISTRY DEPARTMENT ATTN T BRILL NEWARK DE 19711
1	UNIVERSITY OF FLORIDA DEPT OF CHEMISTRY ATTN J WINEFORDNER GAINESVILLE FL 32611
3	GEORGIA INSTITUTE OF TECHNOLOGY SCHOOL OF AEROSPACE ENGINEERING ATTN E PRICE W C STRAHLE B T ZINN ATLANTA GA 30332
1	UNIVERSITY OF ILLINOIS DEPT OF MECH ENG ATTN H KRIER 144MEB 1206 W GREEN ST URBANA IL 61801
1	THE JOHNS HOPKINS UNIV CPIA ATTN T W CHRISTIAN 10630 LITTLE PATUXENT PKWY SUITE 202 COLUMBIA MD 21044-3200
1	UNIVERSITY OF MICHIGAN GAS DYNAMICS LAB ATTN G M FAETH AEROSPACE ENGINEERING BLDG ANN ARBOR MI 48109-2140
1	UNIVERSITY OF MINNESOTA DEPT OF MECHANICAL ENGINEERING ATTN E FLETCHER MINNEAPOLIS MN 55455

NO. OF COPIES	ORGANIZATION
4	PENNSYLVANIA STATE UNIVERSITY DEPT OF MECHANICAL ENGINEERING ATTN K KUO M MICCI S THYNELL V YANG UNIVERSITY PARK PA 16802
2	PRINCETON UNIVERSITY FORRESTAL CAMPUS LIBRARY ATTN K BREZINSKY I GLASSMAN P O BOX 710 PRINCETON NJ 08540
1	PURDUE UNIVERSITY SCHL OF AERONAUTICS & ASTRONAUTICS ATTN J R OSBORN GRISSOM HALL WEST LAFAYETTE IN 47906
1	PURDUE UNIVERSITY DEPARTMENT OF CHEMISTRY ATTN E GRANT WEST LAFAYETTE IN 47906
2	PURDUE UNIVERSITY SCHL OF MECHANICAL ENGNRNG ATTN N M LAURENDEAU S N B MURTHY TSPC CHAFFEE HALL WEST LAFAYETTE IN 47906
1	RENSSELAER POLYTECHNIC INST DEPT OF CHEMICAL ENGINEERING ATTN A FONTIJN TROY NY 12181
1	STANFORD UNIVERSITY DEPT OF MECHANICAL ENGINEERING ATTN R HANSON STANFORD CA 94305
1	UNIVERSITY OF TEXAS DEPT OF CHEMISTRY ATTN W GARDINER AUSTIN TX 78712
1	VA POLYTECH INST AND STATE UNIV ATTN J A SCHETZ BLACKSBURG VA 24061

<u>NO. OF COPIES</u>	<u>ORGANIZATION</u>
1	APPLIED COMBUSTION TECHNOLOGY INC ATTN A M VARNEY P O BOX 607885 ORLANDO FL 32860
2	APPLIED MECHANICS REVIEWS ASME ATTN R E WHITE A B WENZEL 345 E 47TH STREET NEW YORK NY 10017
1	TEXTRON DEFENSE SYSTEMS ATTN A PATRICK 2385 REVERE BEACH PARKWAY EVERETT MA 02149-5900
1	BATTELLE TWSTIAC 505 KING AVENUE COLUMBUS OH 43201-2693
1	COHEN PROFESSIONAL SERVICES ATTN N S COHEN 141 CHANNING STREET REDLANDS CA 92373
1	EXXON RESEARCH & ENG CO ATTN A DEAN ROUTE 22E ANNANDALE NJ 08801
1	GENERAL APPLIED SCIENCE LABS INC 77 RAYNOR AVENUE RONKONKAMA NY 11779-6649
1	GENERAL ELECTRIC ORDNANCE SYSTEMS ATTN J MANDZY 100 PLASTICS AVENUE PITTSFIELD MA 01203
1	GENERAL MOTORS RSCH LABS PHYSICAL CHEMISTRY DEPARTMENT ATTN T SLOANE WARREN MI 48090-9055
2	HERCULES INC ATTN W B WALKUP E A YOUNT P O BOX 210 ROCKET CENTER WV 26726

<u>NO. OF COPIES</u>	<u>ORGANIZATION</u>
1	HERCULES INC ATTN R V CARTWRIGHT 100 HOWARD BLVD KENVIL NJ 07847
1	ALLIANT TECHSYSTEMS INC MARINE SYSTEMS GROUP ATTN D E BRODEN MS MN50 2000 600 2ND STREET NE HOPKINS MN 55343
1	ALLIANT TECHSYSTEMS INC ATTN R E TOMPKINS MN 11 2720 600 SECOND ST NORTH HOPKINS MN 55343
1	IBM CORPORATION RESEARCH DIVISION ATTN A C TAM 5600 COTTLE ROAD SAN JOSE CA 95193
1	IIT RESEARCH INSTITUTE ATTN R F REMALY 10 WEST 35TH STREET CHICAGO IL 60616
1	LOCKHEED MISSILES & SPACE CO ATTN GEORGE LO 3251 HANOVER STREET DEPT 52 35 B204 2 PALO ALTO CA 94304
1	OLIN ORDNANCE ATTN V MCDONALD LIBRARY P O BOX 222 ST MARKS FL 32355-0222
1	PAUL GOUGH ASSOCIATES INC ATTN P S GOUGH 1048 SOUTH STREET PORTSMOUTH NH 03801-5423
1	HUGHES AIRCRAFT COMPANY ATTN T E WARD PO BOX 11337 TUCSON AZ 85734-1337

<u>NO. OF COPIES</u>	<u>ORGANIZATION</u>
1	SCIENCE APPLICATIONS INC ATTN R B EDELMAN 23146 CUMORAH CREST WOODLAND HILLS CA 91364
3	SRI INTERNATIONAL ATTN G SMITH D CROSLEY D GOLDEN 333 RAVENSWOOD AVENUE MENLO PARK CA 94025
1	STEVENS INSTITUTE OF TECH DAVIDSON LABORATORY ATTN R MCALEVY III HOBOKEN NJ 07030
1	SVERDRUP TECHNOLOGY INC LERC GROUP ATTN R J LOCKE MS SVR 2 2001 AEROSPACE PARKWAY BROOK PARK OH 44142
1	SVERDRUP TECHNOLOGY INC ATTN J DEUR 2001 AEROSPACE PARKWAY BROOK PARK OH 44142
3	THIOKOL CORPORATION ELKTON DIVISION ATTN R BIDDLE R WILLER TECH LIB P O BOX 241 ELKTON MD 21921
3	THIOKOL CORPORATION WASATCH DIVISION ATTN S J BENNETT P O BOX 524 BRIGHAM CITY UT 84302
1	UNITED TECHNOLOGIES RSCH CENTER ATTN A C ECKBRETH EAST HARTFORD CT 06108
1	UNITED TECHNOLOGIES CORP CHEMICAL SYSTEMS DIVISION ATTN R R MILLER P O BOX 49028 SAN JOSE CA 95161-9028

<u>NO. OF COPIES</u>	<u>ORGANIZATION</u>
1	UNIVERSAL PROPULSION COMPANY ATTN H J MCSPADEN 25401 NORTH CENTRAL AVENUE PHOENIX AZ 85027-7837
1	VERITAY TECHNOLOGY INC ATTN E B FISHER 4845 MILLERSPORT HIGHWAY EAST AMHERST NY 14051-0305
1	FREEDMAN ASSOCIATES ATTN E FREEDMAN 2411 DIANA ROAD BALTIMORE MD 21209-1525
3	ALLIANT TECHSYSTEMS ATTN C CANDLAND L OSGOOD R BECKER 600 SECOND ST NE HOPKINS MN 55343
1	US ARMY BENET LABORATORY ATTN SAM SOPOK SMCAR CCB B WATERVLIET NY 12189

NO. OF
COPIES ORGANIZATION

ABERDEEN PROVING GROUND

36 DIR USARL
ATTN: AMSRL-WT-P, A HORST
 AMSRL-WT-PC,
 R A FIFER
 G F ADAMS
 W R ANDERSON
 R A BEYER
 S W BUNTE
 C F CHABALOWSKI
 K P MCNEILL-BOONSTOPPEL
 A COHEN
 R CUMPTON
 R DANIEL
 D DEVYNCK
 N F FELL
 B E FORCH
 J M HEIMERL
 A J KOTLAR
 M R MANAA
 W F MCBRATNEY
 K L MCNESBY
 S V MEDLIN
 M S MILLER
 A W MIZIOLEK
 S H MODIANO
 J B MORRIS
 J E NEWBERRY
 S A NEWTON
 R A PESCE-RODRIGUEZ
 B M RICE
 R C SAUSA
 M A SCHROEDER
 J A VANDERHOFF
 M WENSING
 A WHREN
 J M WIDDER
 C WILLIAMSON
 AMSRL-CI-CA, R PATEL

INTENTIONALLY LEFT BLANK.

USER EVALUATION SHEET/CHANGE OF ADDRESS

This Laboratory undertakes a continuing effort to improve the quality of the reports it publishes. Your comments/answers to the items/questions below will aid us in our efforts.

1. ARL Report Number/Author ARL-TR-1232 (Widder) Date of Report November 1996
2. Date Report Received _____
3. Does this report satisfy a need? (Comment on purpose, related project, or other area of interest for which the report will be used.) _____

4. Specifically, how is the report being used? (Information source, design data, procedure, source of ideas, etc.) _____

5. Has the information in this report led to any quantitative savings as far as man-hours or dollars saved, operating costs avoided, or efficiencies achieved, etc? If so, please elaborate. _____

6. General Comments. What do you think should be changed to improve future reports? (Indicate changes to organization, technical content, format, etc.) _____

CURRENT
ADDRESS

Organization

Name

Street or P.O. Box No.

City, State, Zip Code

7. If indicating a Change of Address or Address Correction, please provide the Current or Correct address above and the Old or Incorrect address below.

OLD
ADDRESS

Organization

Name

Street or P.O. Box No.

City, State, Zip Code

(Remove this sheet, fold as indicated, tape closed, and mail.)
(DO NOT STAPLE)


DEPARTMENT OF THE ARMY

OFFICIAL BUSINESS

BUSINESS REPLY MAIL
FIRST CLASS PERMIT NO 0001,APG,MD

POSTAGE WILL BE PAID BY ADDRESSEE

**DIRECTOR
US ARMY RESEARCH LABORATORY
ATTN AMSRL WT PC
ABERDEEN PROVING GROUND MD 21005-5066**

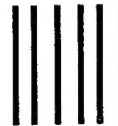

**NO POSTAGE
NECESSARY
IF MAILED
IN THE
UNITED STATES**

